

AD-A237 511



AD

AD-E402 211

Contractor Report ARAED-CR-91008

ANALYSIS OF NITROMETHANE

William R. Herrera
Southwest Research Institute
San Antonio, TX 78228-0510

William O. Seals
Project Engineer

June 1991



US ARMY
ARMAMENT MUNITIONS
& CHEMICAL COMMAND
ARMAMENT RDE CENTER

U.S. ARMY ARMAMENT RESEARCH, DEVELOPMENT AND ENGINEERING CENTER

Armament Engineering Directorate
Picatinny Arsenal, New Jersey

Approved for public release; distribution is unlimited.

REPRODUCED FROM
BEST AVAILABLE COPY

91-03227



The views, opinions, and/or findings contained in this report are those of the author(s) and should not be construed as an official Department of the Army position, policy, or decision, unless so designated by other documentation.

The citation in this report of the names of commercial firms of commercially available products or systems does not constitute official endorsement by or approval of the U.S. Government.

Destroy this report when no longer needed by any method that will prevent disclosure of contents or reconstruction of the document. Do not return to the originator.

REPORT DOCUMENTATION PAGE			Form Approved OMB NO. 0704-0188													
<small>Public reporting burden for this collection of information is estimated to average 1 hour per response, including the time for reviewing instructions, searching existing data sources, gathering and maintaining the data needed, and completing and reviewing the collection of information. Send comments regarding this burden estimate or any other aspect of this collection of information, including suggestions for reducing this burden, to Washington Headquarters Services, Directorate for Information Operations and Reports, 1215 Jefferson Davis Highway, Suite 1204, Arlington, VA 22202-4302, and to the Office of Management and Budget, Paperwork Project (0704-0188), Washington, DC 20503.</small>																
1. AGENCY USE ONLY (Leave blank)	2. REPORT DATE June 1991	3. REPORT TYPE AND DATES COVERED Final, 6/12/90 to 7/31/90														
4. TITLE AND SUBTITLE MIGRAD SYSTEM IMPROVED SAFETY AND PRODUCIBILITY IN THE MANUFACTURING OF PYROTECHNIC MIXTURES		5. FUNDING NUMBERS DAAA21-88-D-0021														
6. AUTHOR(S) William R. Herrera, Southwest Research Institute William O. Seals, ARDEC																
7. PERFORMING ORGANIZATION NAME(S) AND ADDRESS(ES) Southwest Research Institute San Antonio, TX 78228-0510		8. PERFORMING ORGANIZATION REPORT NUMBER ARDEC, AED ATTN: Energetics and Warheads Div ATTN: SMCAR-AES-M Picatinny Arsenal, NJ 07806-5000														
9. SPONSORING/MONITORING AGENCY NAME(S) AND ADDRESS(ES) ARDEC, IMD STINFO Br ATTN: SMCAR-IMI-I Picatinny Arsenal, NJ 07806-5000		10. SPONSORING/MONITORING AGENCY REPORT NUMBER Contractor Report ARAED-CR-91008														
11. SUPPLEMENTARY NOTES																
12a. DISTRIBUTION/AVAILABILITY STATEMENT Approved for public release; distribution is unlimited.		12b. DISTRIBUTION CODE														
13. ABSTRACT (Maximum 200 words) Samples were taken at random from 60 barrels of nitromethane which were stored at Yuma Proving Grounds and sent to Southwest Research Institute for analyses. The nitromethane content was determined for each along with contaminate analyses for nitroethane, 2-nitropropane, water, and metals. The results indicated the following:																
<table border="1"> <thead> <tr> <th>Component</th> <th>Percent</th> </tr> </thead> <tbody> <tr> <td>Nitromethane</td> <td>98 or greater</td> </tr> <tr> <td>Nitroethane</td> <td>1.1 or less</td> </tr> <tr> <td>2-Nitropropane</td> <td>0.2 or less</td> </tr> <tr> <td>Water</td> <td>0.11 or less</td> </tr> <tr> <td>Metals</td> <td>Trace (less than 10 ppb)</td> </tr> </tbody> </table>					Component	Percent	Nitromethane	98 or greater	Nitroethane	1.1 or less	2-Nitropropane	0.2 or less	Water	0.11 or less	Metals	Trace (less than 10 ppb)
Component	Percent															
Nitromethane	98 or greater															
Nitroethane	1.1 or less															
2-Nitropropane	0.2 or less															
Water	0.11 or less															
Metals	Trace (less than 10 ppb)															
14. SUBJECT TERMS Nitromethane; Nitroethane; FTIR; Karl Fischer; Volatiles; Semivolatiles; Mallinckrodt gas chromatography; Atomic absorption			15. NUMBER OF PAGES 61													
			16. PRICE CODE													
17. SECURITY CLASSIFICATION OF REPORT UNCLASSIFIED	18. SECURITY CLASSIFICATION OF THIS PAGE UNCLASSIFIED	19. SECURITY CLASSIFICATION OF ABSTRACT UNCLASSIFIED	20. LIMITATION OF ABSTRACT SAR													

CONTENTS

	Page
Introduction	1
Procedures and Results	2
Discussion	6
Appendixes	
A FTIR Comparison Spectra	7
B FTIR Spectra of Samples	13
C Volatile Analyses	25
Distribution List	47

TABLES

1 FTIR analysis for percent nitromethane	2
2 GC/MS analysis for volatile contaminants	3
3 GC/MS analysis for semivolatile contaminants in selected nitromethane samples	3
4 Karl Fischer analysis of nitromethane samples	4
5 Metals analysis of nitromethane samples	4
6 Computation of various analysis of nitromethane samples	5

Accession For	
NTIS GRANT	<input checked="" type="checkbox"/>
DTIC TAB	<input type="checkbox"/>
Unannounced	<input type="checkbox"/>
Justification	
By	
Distribution/	
Availability Codes	
Dist	Avail and/or Special
A-1	

1. INTRODUCTION

Originally, the Government had purchased and stored 60 barrels of nitromethane at the Yuma Proving Grounds for use in the TEXS system. This system required that the materials be housed in the all-purpose storage area. DOD has classified nitromethane as a Class 4, mass detonating material and not suitable for storage in the all-purpose area. This restriction resulted in the rejection of nitromethane as a candidate for the TEXS system.

These barrels have been exposed under varying degrees of environmental conditions and changes. Contamination and degradation to the nitromethane may have occurred during the storage period. An analysis is required to determine whether the material has deviated dramatically from the specifications furnished by the supplier. This analysis will determine whether the stored nitromethane is suitable for use in another application.

Samples of nitromethane were taken at random from the 60 barrels of stored nitromethane. These samples were taken from the top, middle, and bottom to ascertain the uniformity of nitromethane composition. Replicate samples of each were analyzed in the following manner for these components and contaminants that may be present:

Component	Analytical Method
Nitromethane	FTIR
Nitroethane	GC/Mass Spectroscopy
2-Nitropropane	GC/Mass Spectroscopy
Water	Karl Fischer
Metals	AA Graphite Furnace

FTIR - Fourier transform infrared

GC - gas chromatography

AA - atomic absorption

2. PROCEDURES AND RESULTS

Analysis of nitromethane was performed by FTIR. SwRI received a 99%+ sample of nitromethane that was certified by Angus Chemical. An FTIR analysis was performed and the results compared to the Chem Sources 98%+ sample used during the program. The results confirmed the purity level of the Chem Sources sample. The results are given in Table 1 and raw data are found in Appendix A.

TABLE 1. FTIR ANALYSIS FOR PERCENT NITROMETHANE	
Sample	% Nitromethane
Y1-B	99.6
Y1-M	98.3
Y1-T	100.0
Y1-MM	99.3
A1-B	99.3
A1-M	99.9
A1-MM	99.0
N1-B	100.2
N1-M	99.5
N1-T	99.5
Standard	98%

B Bottom
M Middle
T Top
MM Middle mixed

The FTIR spectra of the neat nitromethane samples were compared to that of the Mallinckrodt sample and found to be virtually identical. Computer subtraction of the Mallinckrodt sample from the other samples failed to show any gross contaminant. It was concluded that any impurities present were only in low concentration.

Weighed quantities of each of the samples were dissolved in carbon tetrachloride to make 25 mL total volume. A solution prepared similarly from Chem Service nitromethane (98.0%) was utilized as a standard. This standard was used due to the higher degree of purity of the material. The nitromethane concentration of the samples was determined by comparing absorption intensities at 657 cm^{-1} in a sealed liquid cell. The 657 cm^{-1} absorption is specific for nitromethane and is unaffected by nitroethane and 2-nitropropane. The precision was estimated to be ca. ± 1 nitromethane.

An overlay, offset comparison spectra of Y1-M and the Chem Service 98% standard is given in Appendix A. The two spectra are identical except for a few insignificant peaks in the 800-830 range. The spectra for all the samples are included in Appendix B.

Analysis of volatile compounds other than nitromethane was performed using GC/MS. The results for nitroethane and 2-nitropropane are given in Table 2 and raw data are found in Appendix C.

TABLE 2. GC/MS ANALYSIS FOR VOLATILE CONTAMINANTS		
Sample	% Composition	
	Nitroethane	2-Nitropropane
Y1-B	0.9	~0.1
Y1-M	0.9	~0.1
Y1-T	0.8	~0.1
Y1-MM	0.9	~0.1
A1-B	~0.1	~0.1
A1-M	~0.1	~0.1
A1-MM	~0.4	0.2
N1-T	~0.1	~0.1
N1-M	~0.1	~0.1
N1-B	1.1	~0.1

See legend Table 1.

Analysis for semivolatile compounds was performed using GC/MS. Two samples were selected for analysis. The results are given in Table 3.

TABLE 3. GC/MS ANALYSIS FOR SEMIVOLATILE CONTAMINANTS IN SELECTED NITROMETHANE SAMPLES		
Component	Sample Y1-M ppm	Sample N1-M ppm
Paraldehyde	70.0	38.0
2-Pentanone-4-methyl	42.0	0.0
Total other unknowns	13.4	2.9
% Total Semivolatiles	0.01%	0.004%

See legend Table 1.

No single semivolatile contaminant was found to have a concentration $>0.007\%$.

Karl Fischer analyses were performed to assay the water content of each of the samples. The results can be compared to standard 95% nitromethane sample from Mallinckrodt. The results are given in Table 4.

TABLE 4. KARL FISCHER ANALYSIS OF NITROMETHANE SAMPLES		
Sample	ppm H ₂ O	%
Standard	1,370	0.14
Y1-T	629	0.06
Y1-MM	697	0.07
Y1-M	1077	0.11
Y1-B	693	0.07
N1-T	705	0.07
N1-M	506	0.05
N1-B	432	0.04
A1-M	549	0.05
A1-MM	589	0.06
A1-B	590	0.06

See legend Table 1.

The average percent water of the samples was 0.06.

Metals contamination analysis was performed on a Perkin-Elmer 5000 graphite furnace atomic absorption spectrophotometer. Calibration standards were made up in methanol and samples were diluted 1:1 in methanol before analysis. A methanol blank was run with each sample set. No solvent contamination was observed.

The results are presented in Table 5.

TABLE 5. METALS ANALYSIS OF NITROMETHANE SAMPLES			
Sample	ppb		
	Cu	Fe	Ni
Standard	<1.0	<2.0	<5.0
Y1-T	8.2	<4.0	<10.0
Y1-MM	9.0	<4.0	<10.0
Y1-M	9.0	<4.0	<10.0

TABLE 5. METALS ANALYSIS OF NITROMETHANE SAMPLES			
Sample -	ppb		
	Cu	Fe	Ni
Y1-B	6.6	<4.0	<10.0
N1-T	6.6	<4.0	<10.0
N1-M	6.6	<4.0	<10.0
N1-B	5.8	<4.0	<10.0
A1-MM	9.8	<4.0	<10.0
A1-M	8.0	<4.0	<10.0
A1-B	9.8	<4.0	<10.0

See legend Table 1.

A compendium of all the analytical results are given in Table 6.

TABLE 6. COMPUTATION OF VARIOUS ANALYSIS ON NITROMETHANE SAMPLES						
Sample	% Composition					Metal Contaminants (ppb)
	Nitromethane	H ₂ O	Nitroethane	2-Nitropropane	Semivolatiles	
Y1-B	99.6	0.07	0.9	~0.1		<10
Y1-M	98.3	0.11	0.9	~0.1	0.01	<10
Y1-T	100.0	0.06	0.8	~0.1		<10
Y1-MM	99.3	0.07	0.9	~0.1		<10
A1-B	99.3	0.06	~0.1	~0.1		<10
A1-M	99.9	0.05	~0.1	~0.1		<10
A1-MM	99.0	0.06	0.4	0.2		<10
N1-B	100.2	0.04	1.1	~0.1		<10
N1-M	99.5	0.05	~0.1	~0.1	<0.01	<10
N1-T	99.5	0.07	~0.1	~0.1		<10

See legend Table 1.

3. DISCUSSION

A high degree of purity of nitromethane was found for all the samples. The FTIR analysis showed that all the samples had a concentration of 98% or greater nitromethane. This finding is based on the assumption that the Chem Service nitromethane standard used for comparison is of 98% purity. Other peaks in the FTIR spectra were of insufficient prominence to quantify, so alternate analytical techniques were utilized to identify the other compounds present.

Analysis of other volatile and semivolatiles, especially nitroethane and 2-nitropropane, was accomplished utilizing GC/MS. It was ascertained that all the samples contained 1.1% or less of nitroethane and 0.2% or less of 2-nitropropane. All other peaks in the volatile analysis were significantly less prominent indicating much lower concentrations of any other volatile contaminants.

The semivolatile analysis yielded only trace quantities of contamination of the nitromethane. Only two samples were selected for analysis due to the similarity of the samples preceding analysis and the small percentage of unaccounted for components. Quantification of contaminants was accomplished by subtracting the effects of the nitromethane. It is possible that a higher percentage of paraldehyde is present due to the elution interference of the solvent, but the results of the prior analysis make this unlikely.


Karl Fisher water analysis yielded results of 0.11% water or less. This is consistent with the results of the metals analysis, which showed only trace amounts of metals contamination. For nitromethane to be corrosive to a steel container, the H_2O concentration must exceed 0.2%. The H_2O concentration was below this limit; hence, there was no or little interaction with the containers.

Sampling procedures did not accompany the samples. The method used to sample the barrels could have an effect on the analytical findings. It will be important in the future to ensure that a standard method of sampling is used, preferably a technique that yields a representative sample from each of the layers within a single drum of nitromethane.


It will also be necessary to establish a chain-of-custody record to ensure proper handling of the samples from the field to the analyst. Such a record could help to facilitate understanding any gross contamination that appears, that are not accounted for in future analysis.


Two possible explanations exist for the discrepancy between the original reported purity of the nitromethane, and the purity found in this program. First, during telephone conversations with Angus Chemical, the original shipment of nitromethane was certified at 96%+. This designation commonly results in product which is percentage points higher. Second, and most probable, if a point sampling procedure was utilized, as opposed to a cross-sectional sampling procedure, then a higher degree of homogeneity would result. Point sampling techniques miss any differences in product due to stratification.

APPENDIX A
FTIR COMPARISON SPECTRA

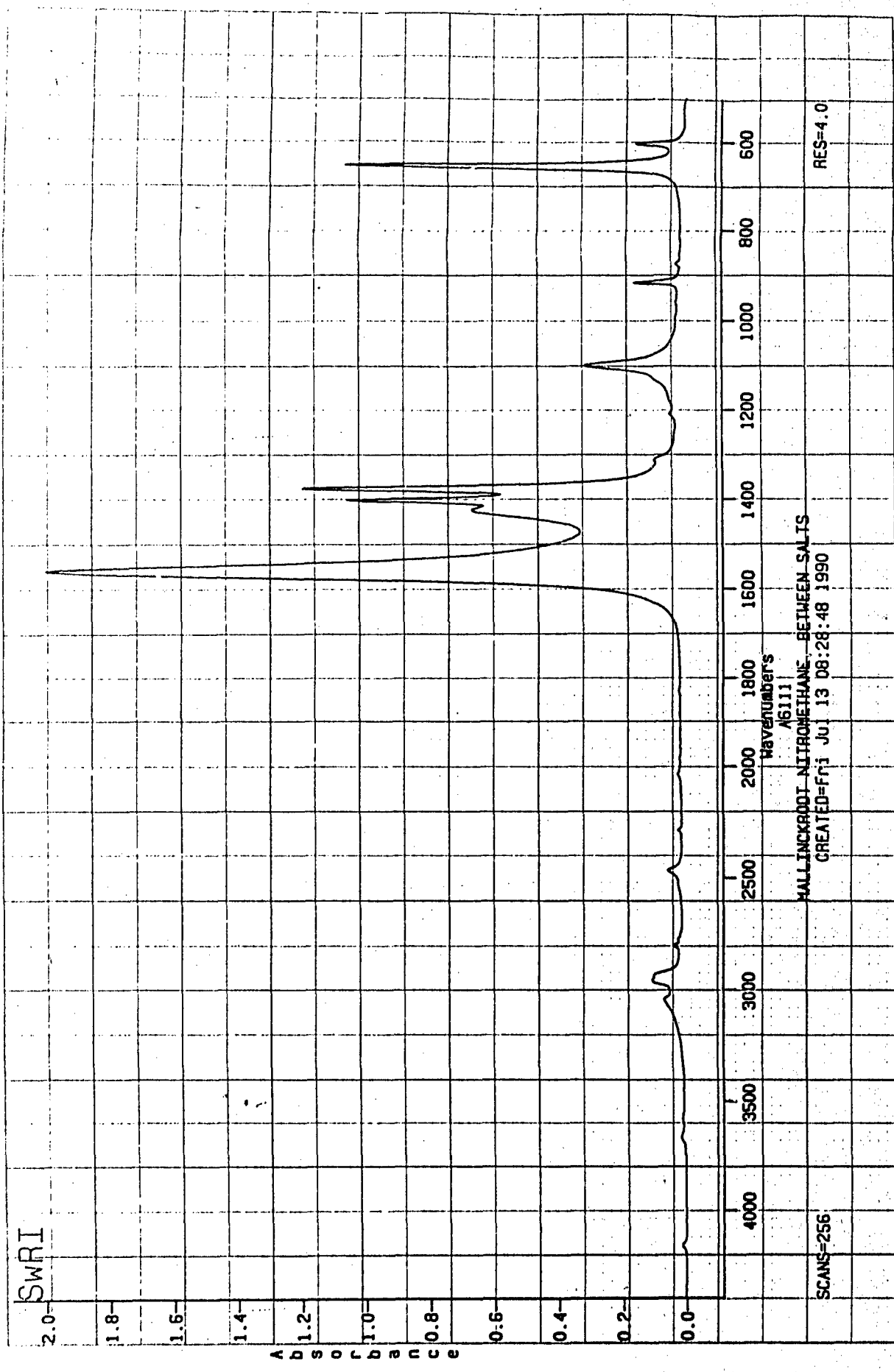

NO. XY 1101-SP5
GALAXY 1101-SP5
PRINTED IN U.S.A.

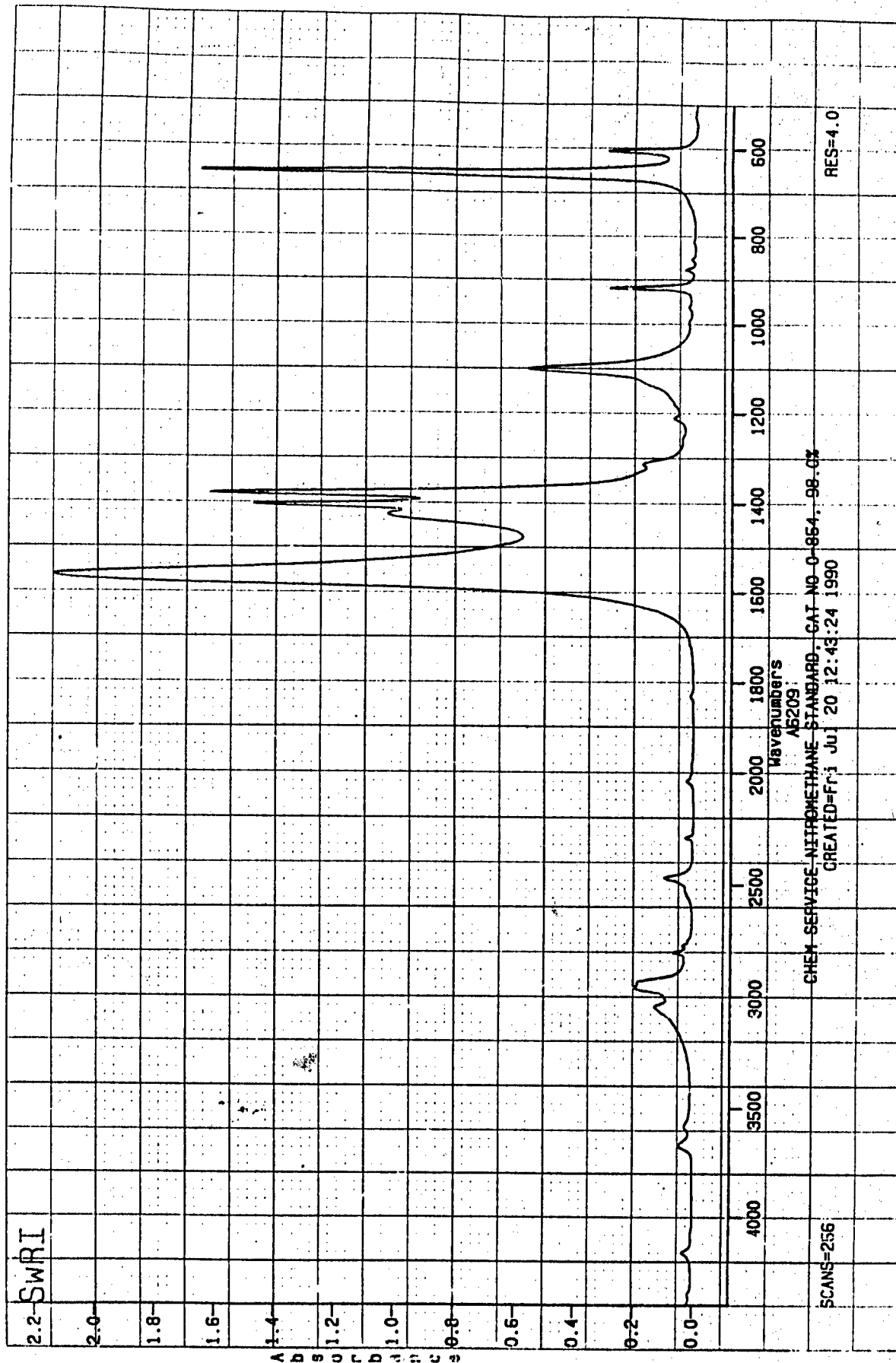


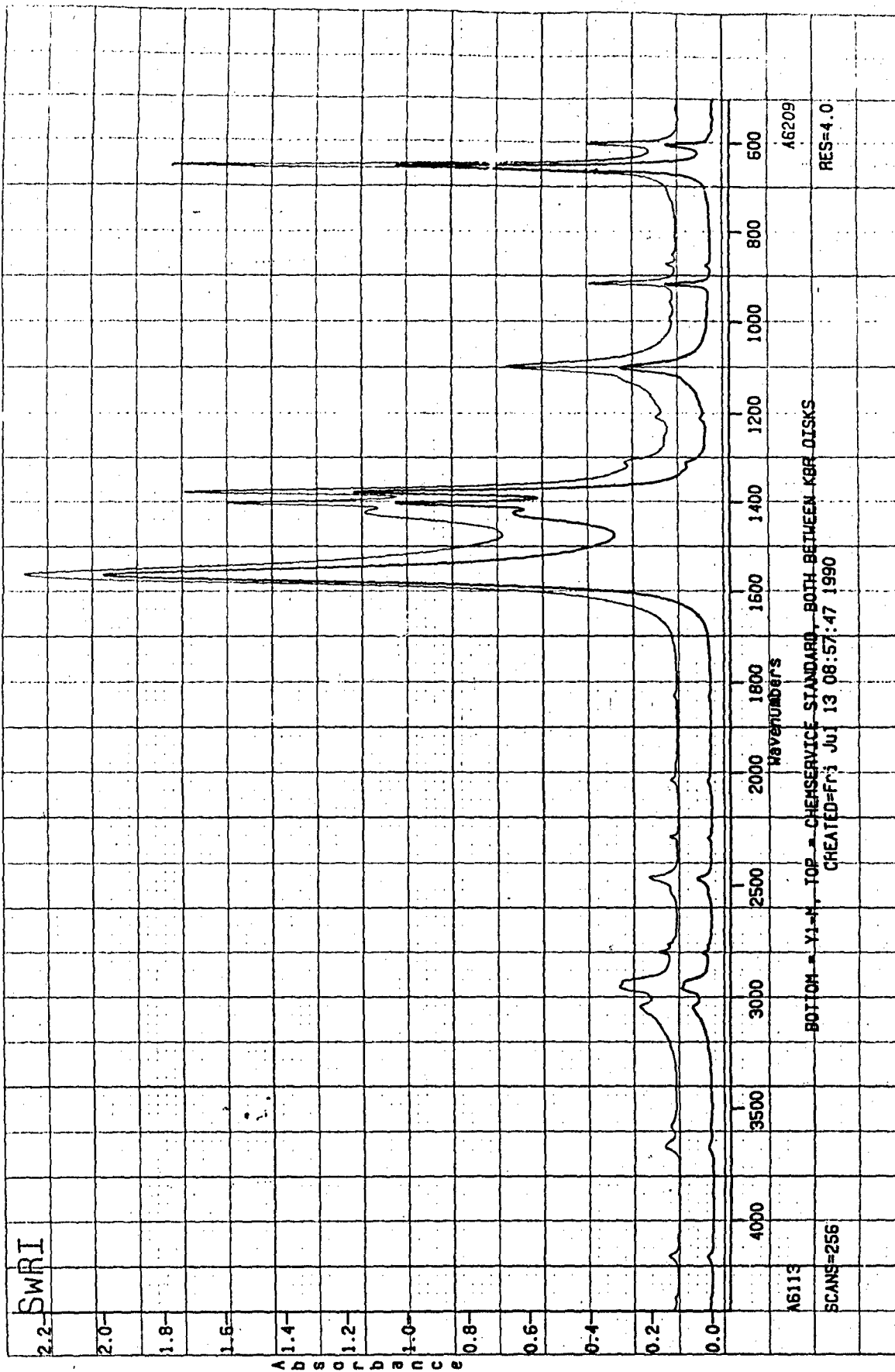

NO. XY 1101-SPS


GRAPHIC SYSTEMS CORPORATION
BURLINGAME, CALIF. 94010
PRINTED IN U.S.A.

APPENDIX B
FTIR SPECTRA OF SAMPLES

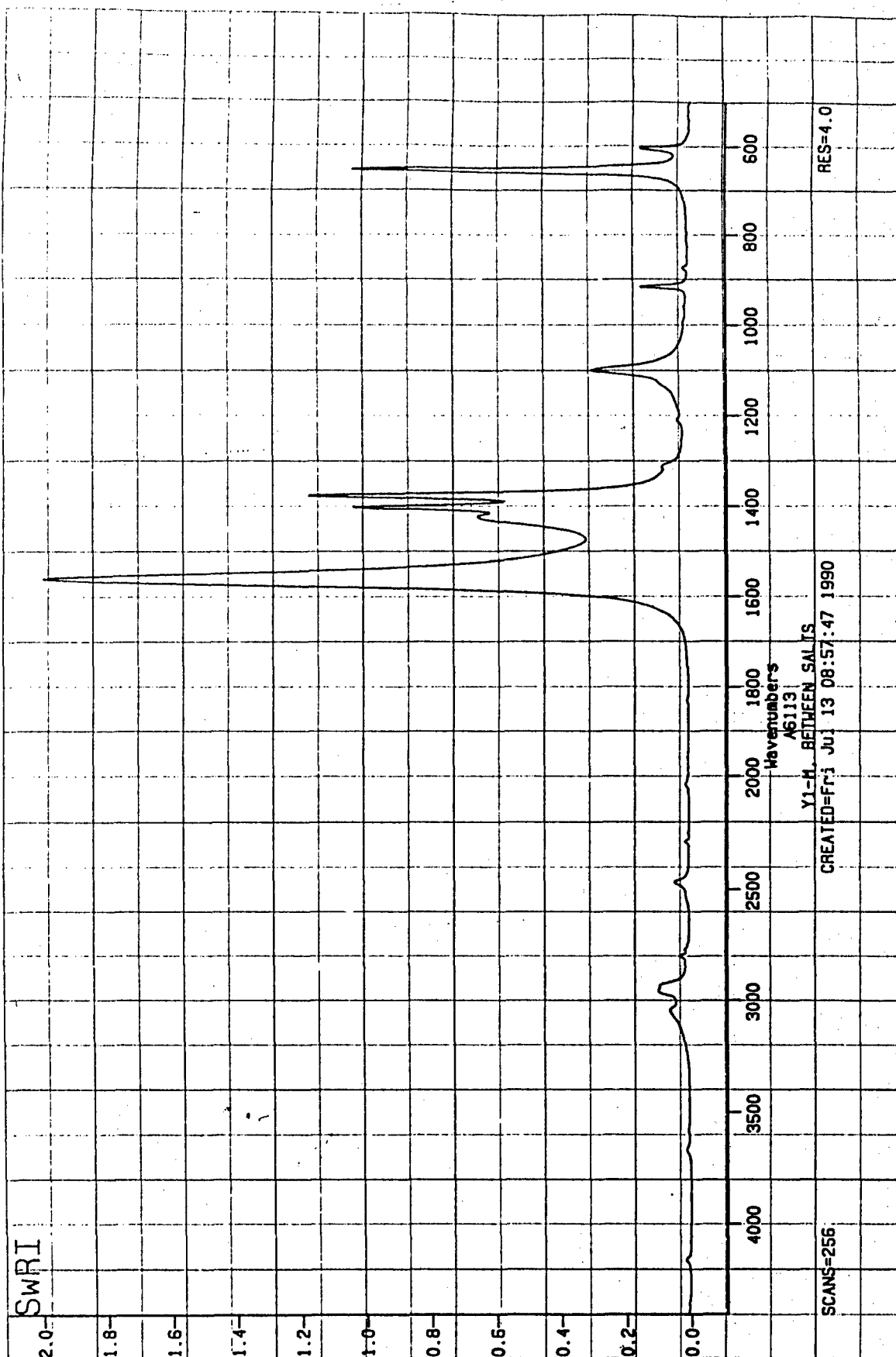


NO. XY 1101-SP5

PERKINELMER
SPECTRA 1000-10000
MODEL 1000-10000
MADE IN U.S.A.

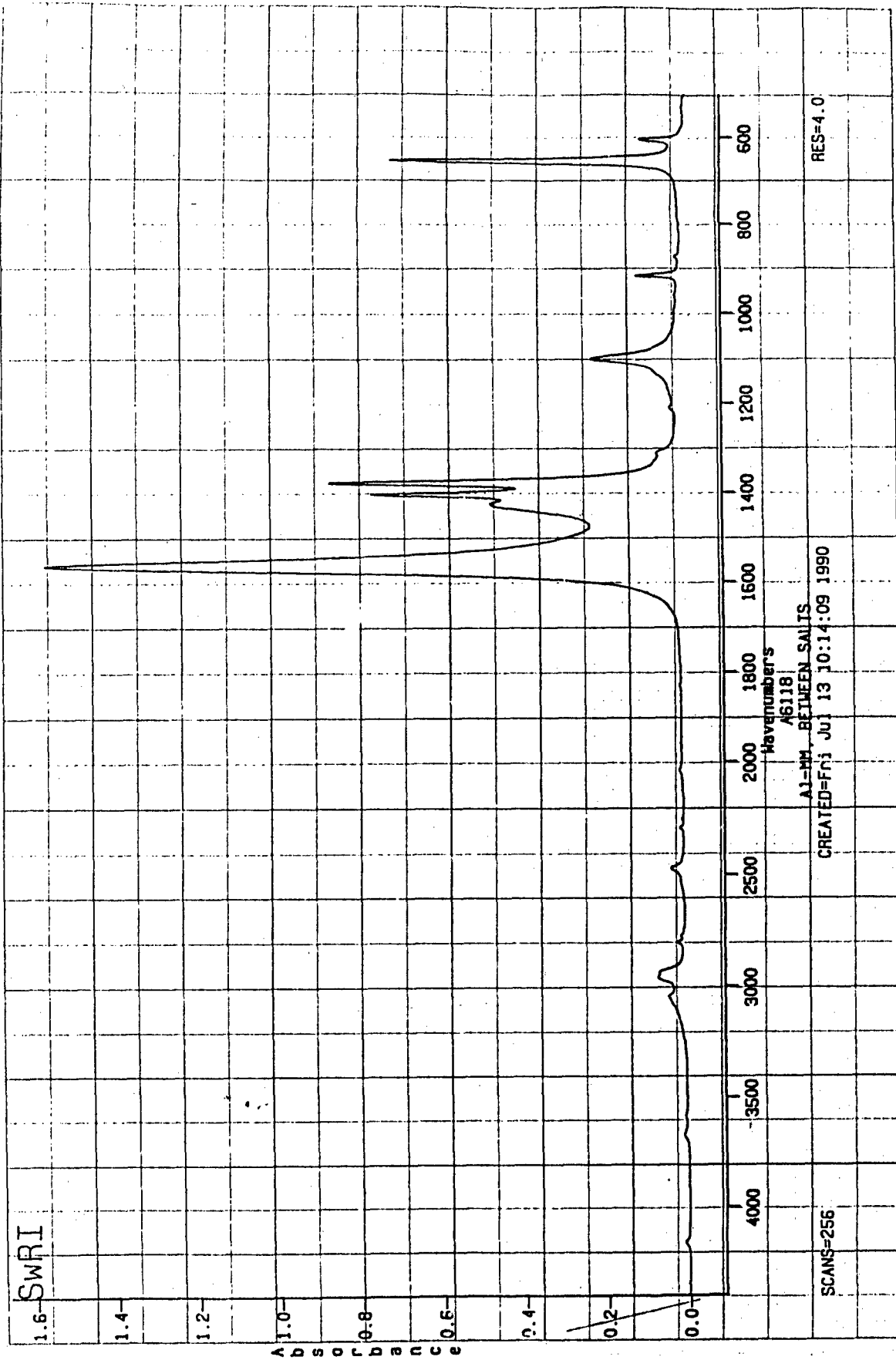

NO XY 1101-SP5

SAFETY PRECAUTIONS
DO NOT TOUCH
THIS EQUIPMENT
UNLESS YOU ARE
PROPERLY TRAINED
AND AUTHORIZED
TO DO SO

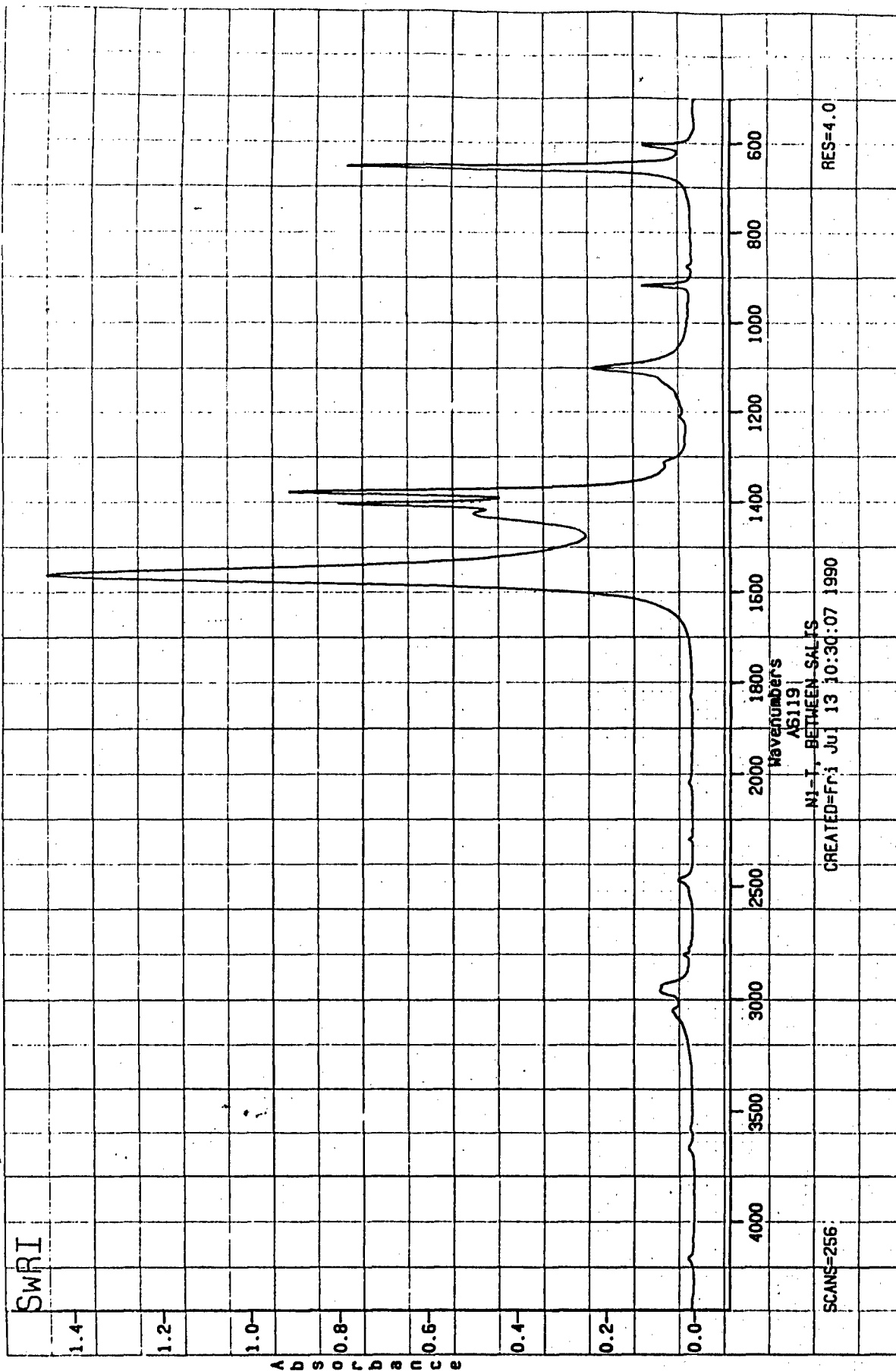

SCANS=256


CREATED=Fri Jul 13 10:30:07 1990

N1-T, BETWEEN SALTS

AS119

Wavenumbers

RES=4.0

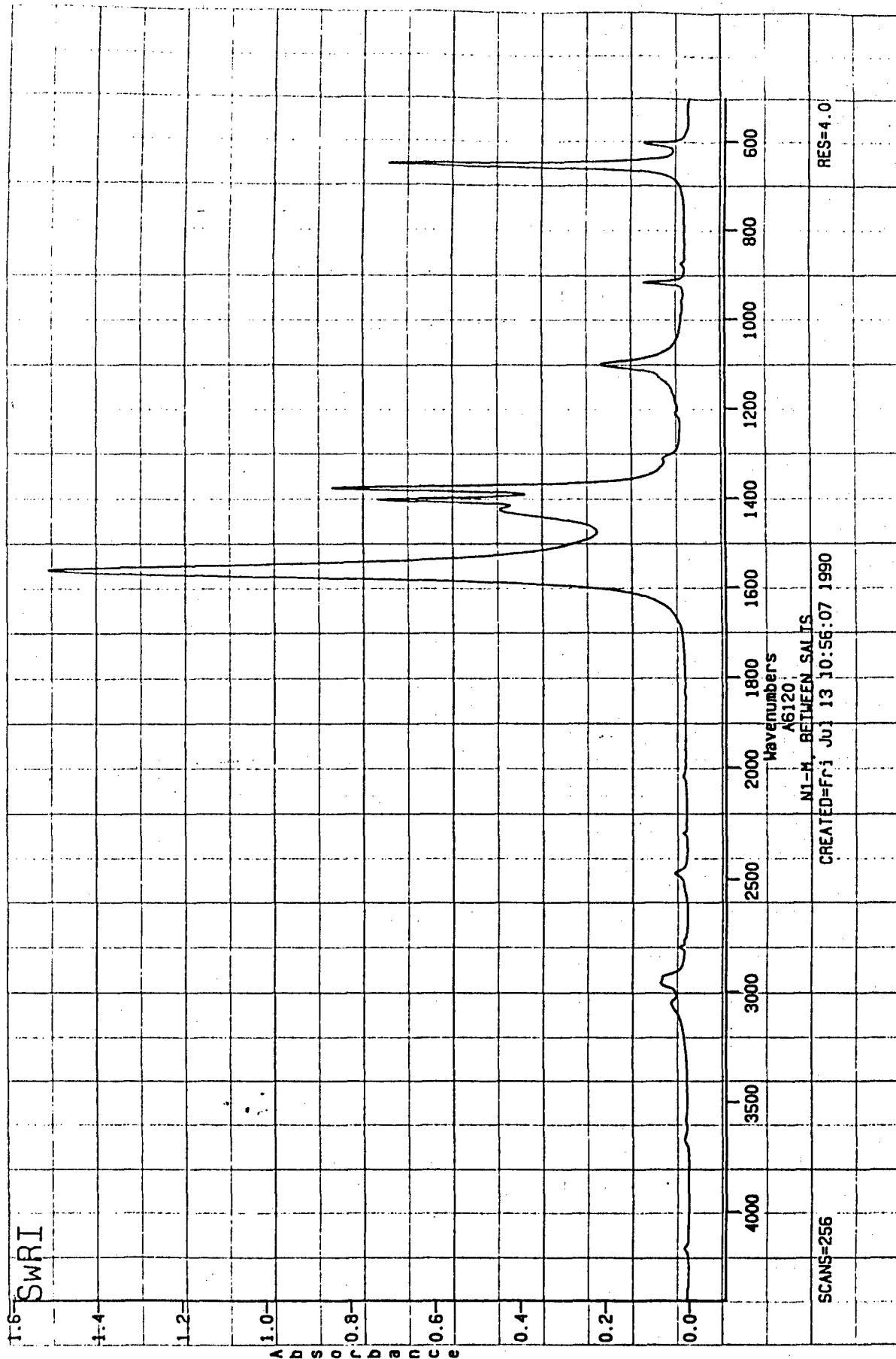


SCANS=256


NI-H, BETWEEN SALTS
A6120
CREATED=Fri Jul 13 10:56:07 1990


RES=4.0

NO XY 1101-SPS

• CALIBRATION CHECK •
LUNDAUL INSTRUMENTS CORPORATION
BURLINGTON, MA 01803
PRINTED IN U.S.A.

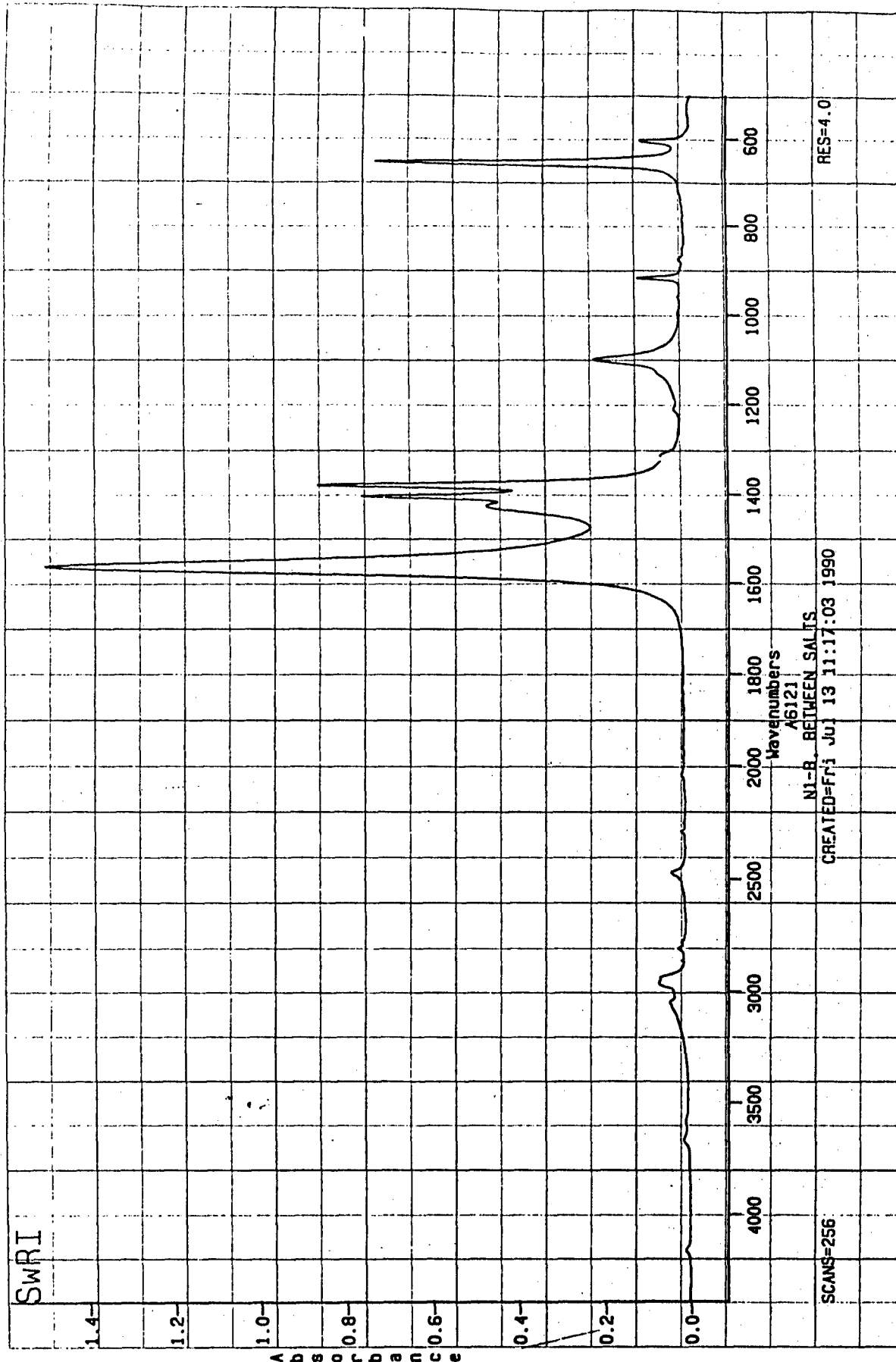

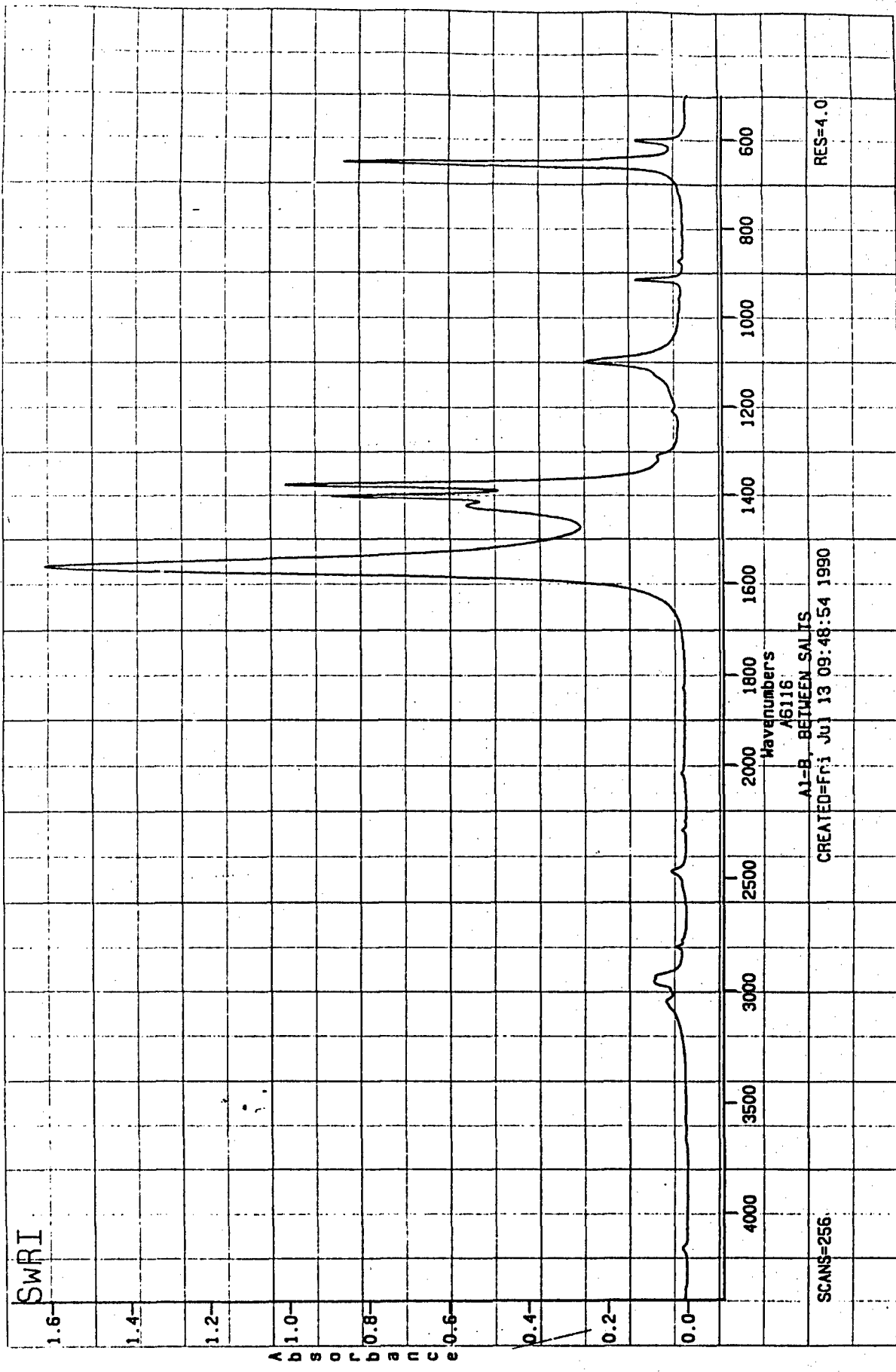


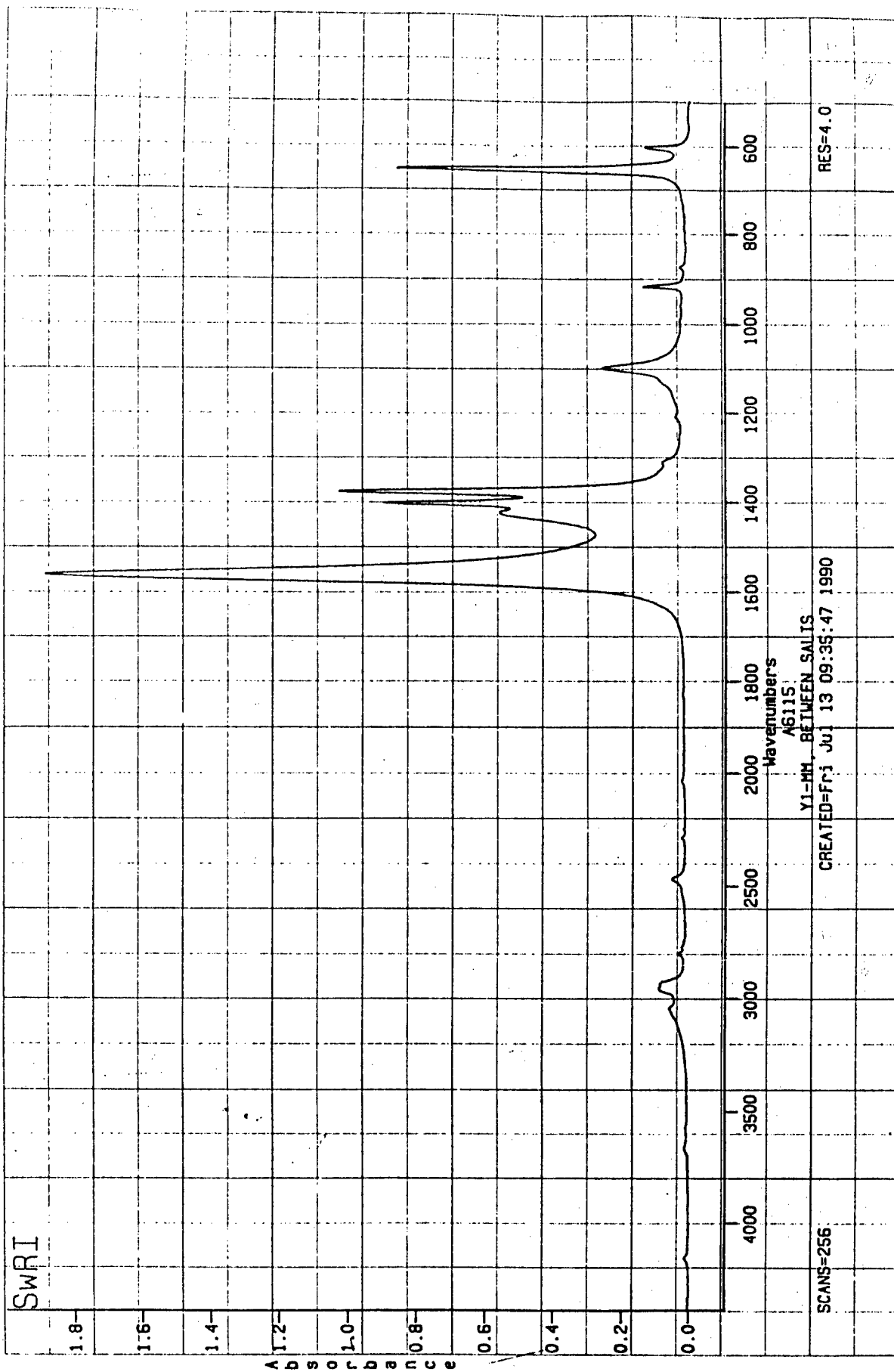

NO XY 1101-SP5

PRINTED IN U.S.A.
SPECTRA LABS & EQUIPMENT, INC.
• 1101-1101-SP5 •

1.8-SWRI

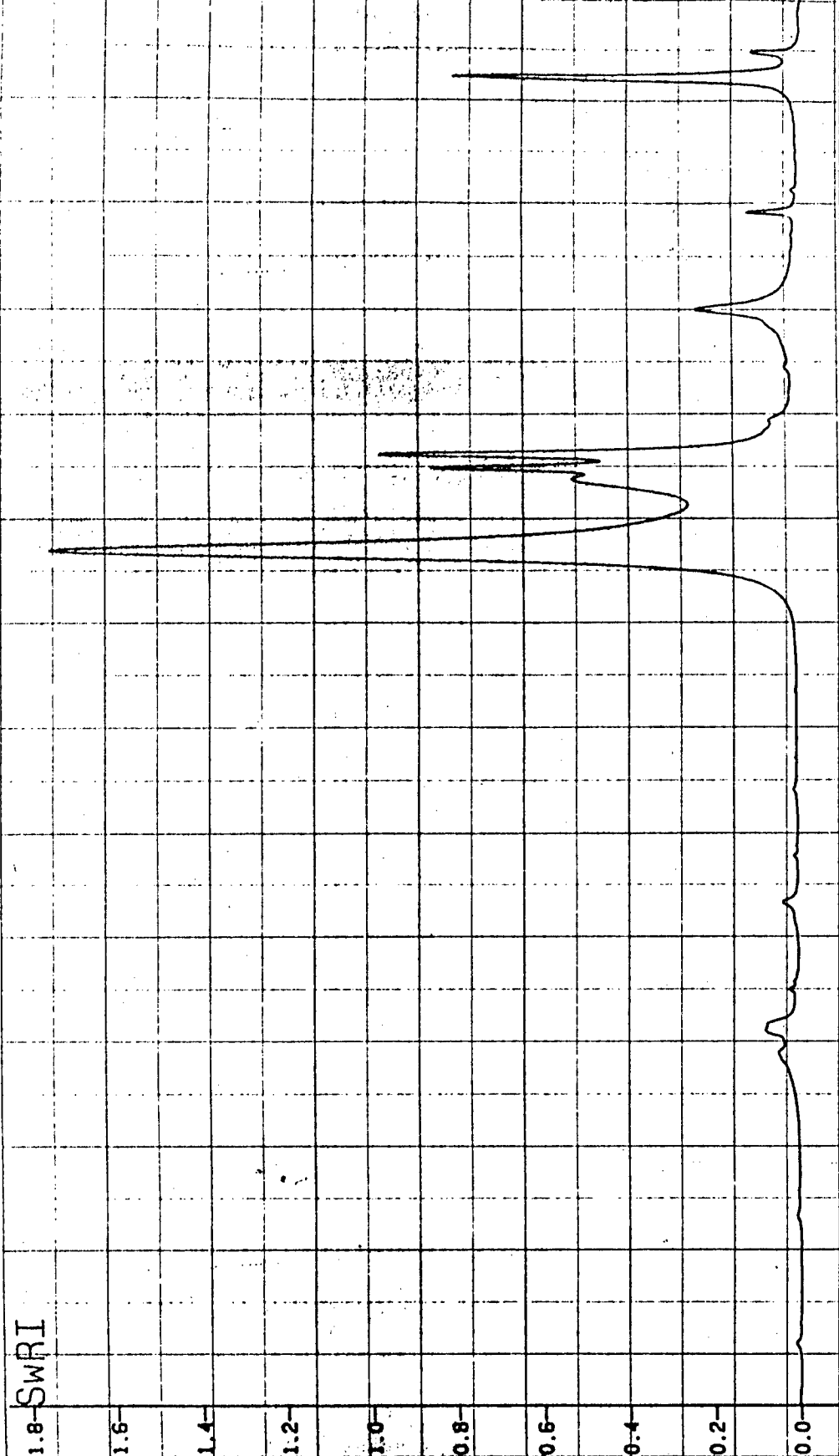

Wavenumbers


AS114

Y1-I, BETWEEN SALIS

CREATED=Fri Jul 13 09:12:00 1990

SCANS=256

RES=4.0

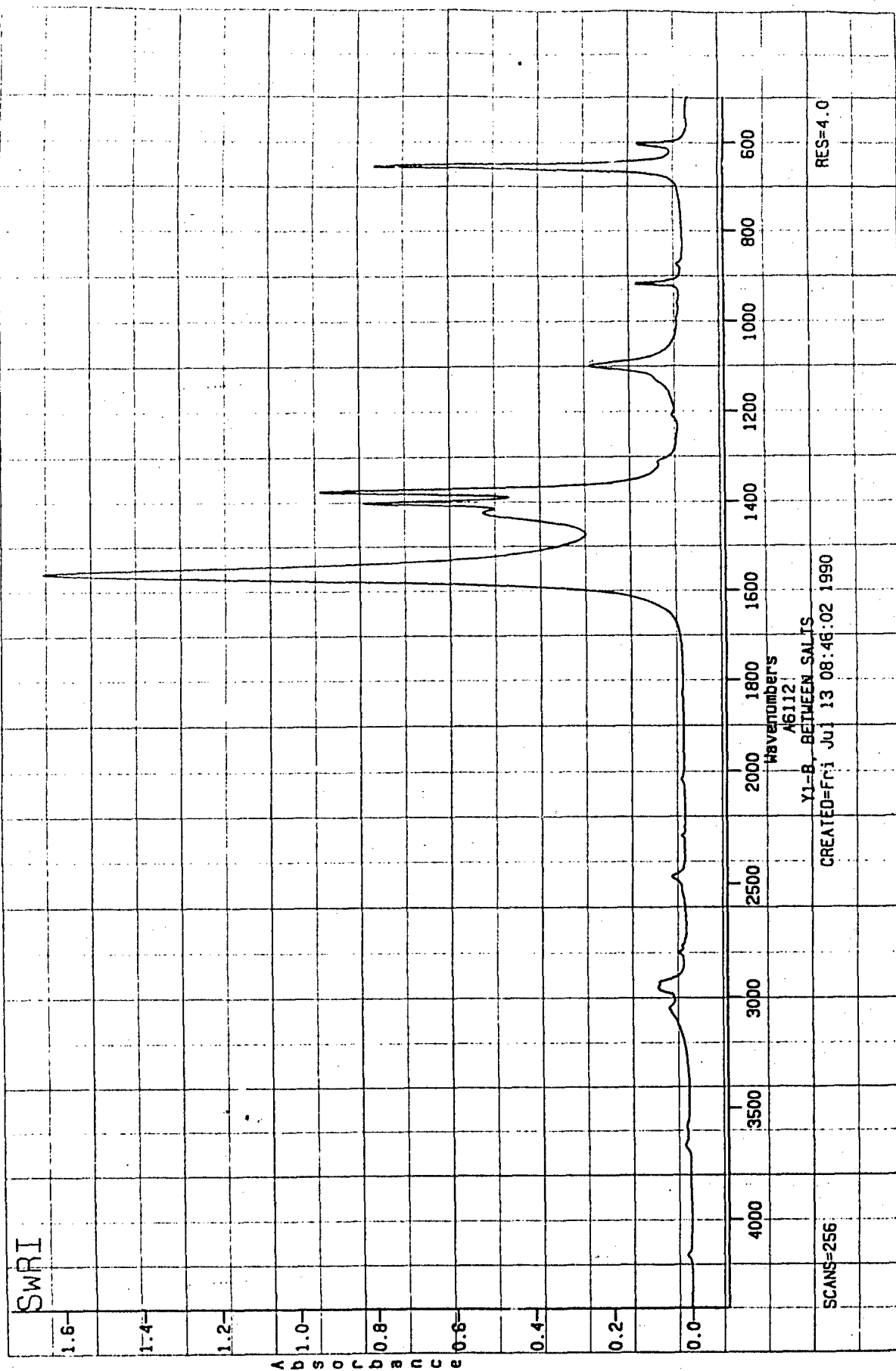

NO XY 1101-SP5

APPENDIX C
VOLATILE ANALYSES


RIC
 05/21/98 23:01:03
 DATA: 00521051 #438
 CALL: 005150F1 #3
 SAMPLE: NITROMETHANE/2-NITROPROPANE 500MG STD: 1UL PURGE
 CONDS.: 300DEG 4MIN/4DEG/MIN/FINAL 180 DEG
 RANGE: G 1, 533 LABEL: N 0, 4.0 QUAN: A 0, 1.0 J 0 BASE: U 20, 3
 434

1724410.

100.07

27


RIC

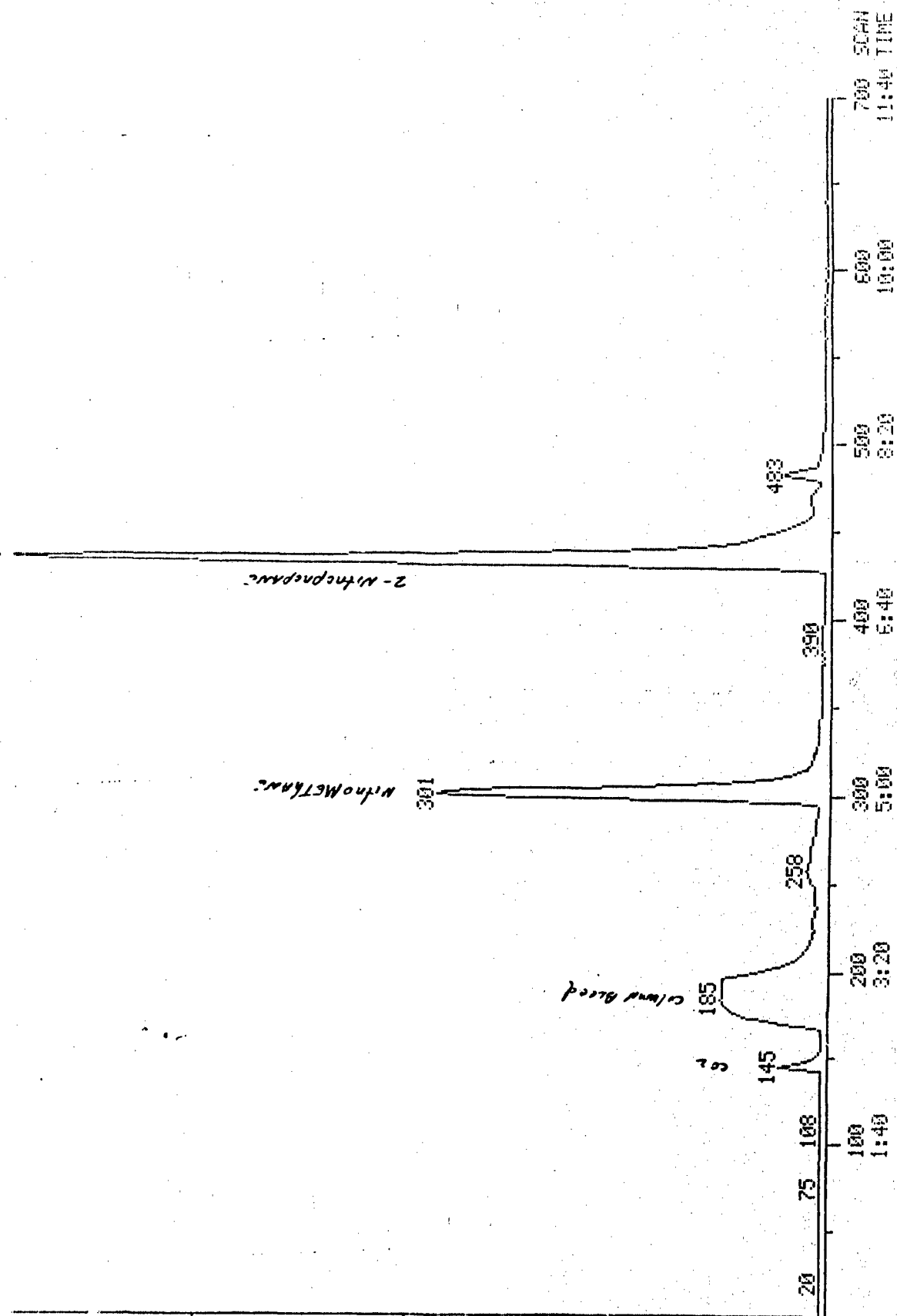

700 SCAN
11:46 TIME

RIC
 06/21/90 20:01:00
 DATA: 00621051 #438
 CALI: 006150F1 #3
 SAMPLE: NITROMETHANE/2-NITROPROPANE 5000NG STD: 1UL PURGE
 COND.: 300DEG 4MIN/40DEG/MIN/FINAL 180 DEG
 RANGE: G 1, 533 LABEL: N 0, 4.0 QUAN: A 0, 1.0 J 0 BASE: U 20, 3
 434

1724410.

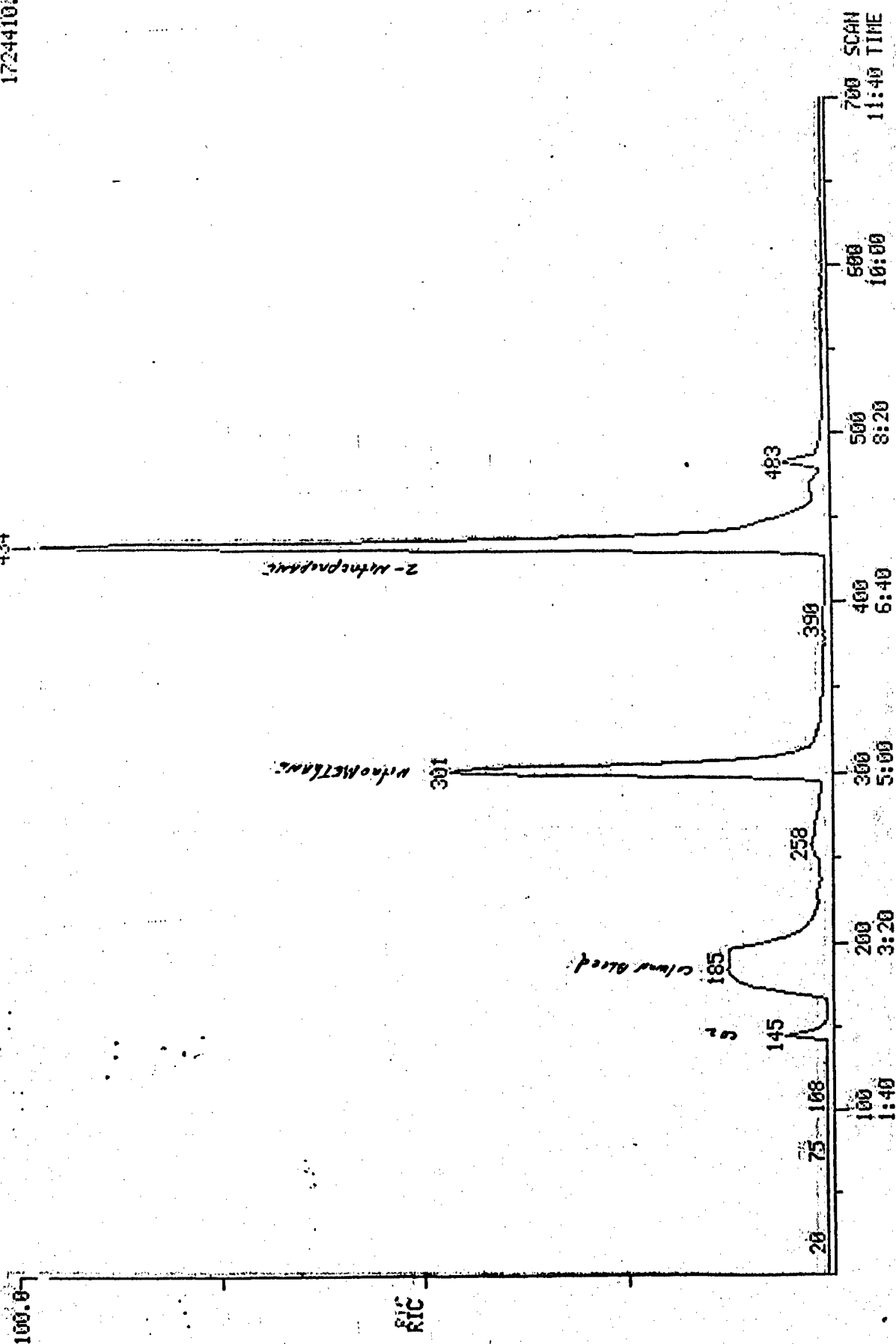

RIC

05/21/90 21:34:00

SAMPLE: UOA METHOD BLANK: 1UL MECH PURGE

CONDOS.: 30DEG 4MIN/40DEG/MIN/FINAL 130 DEG


RANGE: G 1, 533 LABEL: N 0, 4.0 CURN: A 0, 1.0 J 0 BASE: U 20, 3

DATA: C0621081 #452

CALI: C06150F1 #3

SCANS 1 TO 700

100.0

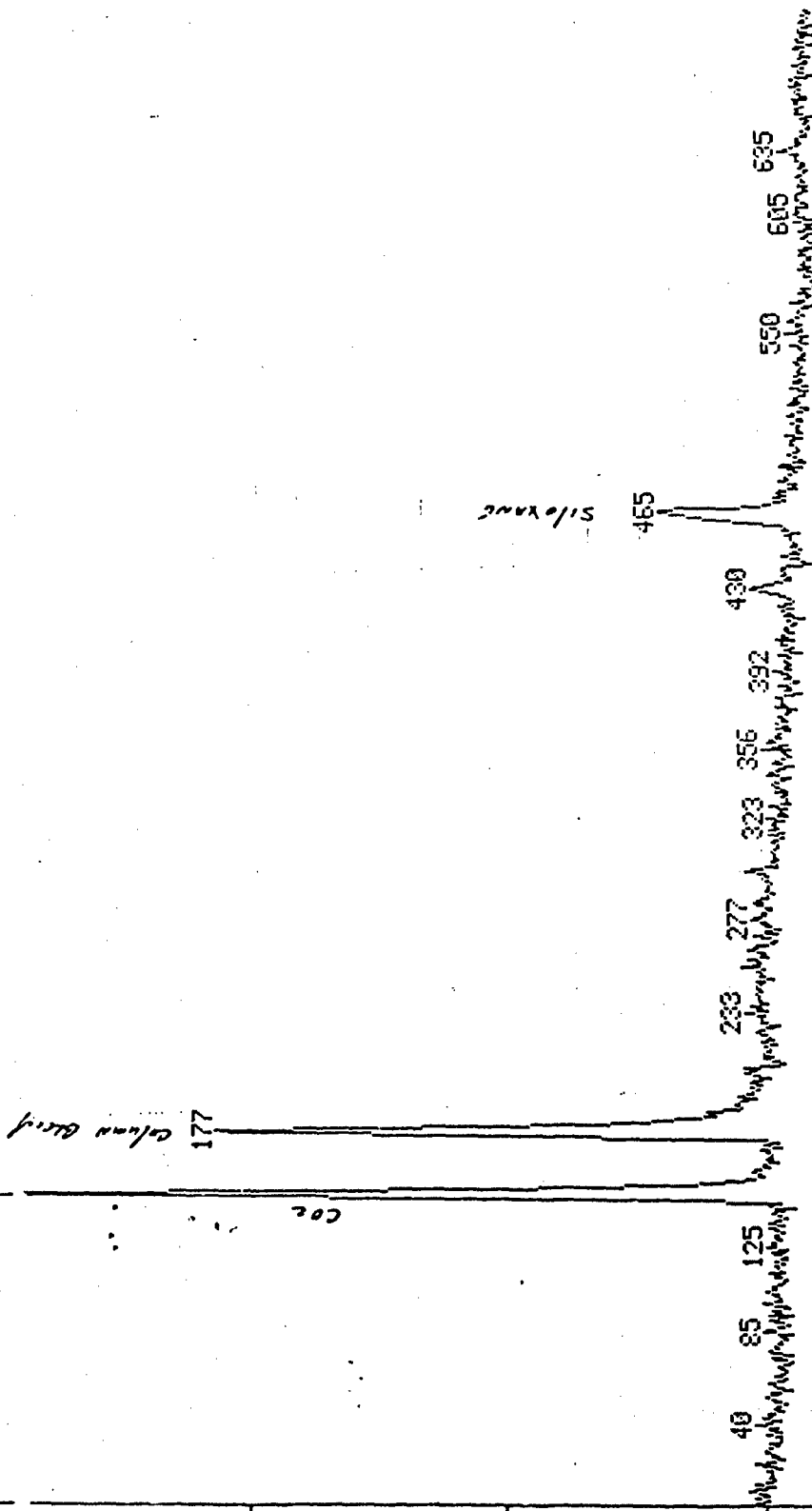

LITERARY SEARCH


05/21/90 21:34:00 + 7:44

SAMPLE: UOA METHOD BLANK: 1UL MEDH PURGE

CONDS.: 30DEG 4MIN/40EG/MIN/FINAL 180 DEG

452 TO # 466 SUMMED - # 470 TO # 475 - # 432 TO # 443 X1.00

DATA: 00521081 # 454
CALI: 005150F1 # 3
BASE M/2: 207
RIC: 22915.

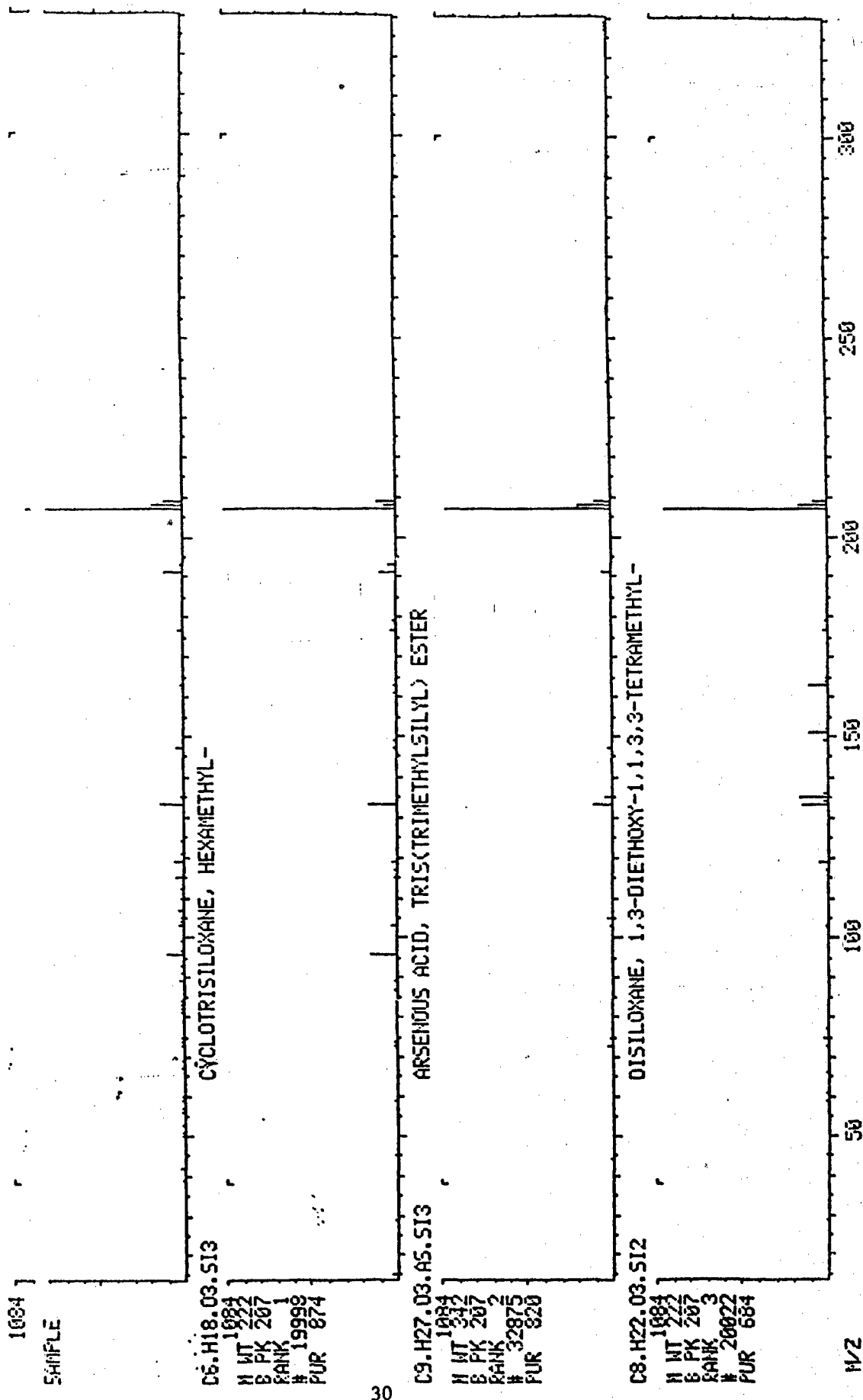

LIESARD, SENFCH

08/21/90 21:34:00 + 7:44

SAMPLE: VOA METHOD BLANK: 1UL MEOH-PURGE

CONDS.: 38DEG 4MIN/40DEG/MIN/FINAL 130 DEG

462 TO # 466 SUMMED - # 470 TO # 475 - # 482 TO # 443 X1.00

DATA: 00621031 # 48+
CALI: 00615061 # 3
R10: 22315.

1034

SAMPLE

08.H18.03.S13

CYCLOTRISILOXANE, HEXAMETHYL-

1984
M WT 223
B PK 207
RANK 1
19998
PUR 874

31

09.H27.03.AS.S13

ARSENIOUS ACID, TRIS(TRIMETHYLSILYL) ESTER

1984
M WT 342
B PK 207
RANK 2
32875
PUR 820

08.H22.03.S12

DISILOXANE, 1,3-DIETHOXY-1,1,3,3-TETRAMETHYL-

1984
M WT 222
B PK 207
RANK 3
20022
PUR 684

H/2

50 100 150 200 250 300

RIC
06/21/98 22:21:00 DATA: C0621001 #427 SCANS 1 TO 700
CALL: C06150F1 #3


SAMPLE: B.H.:A1-B (SI):1UL (50UL SAMPLE/10ML MEDIUM)

CONDUS.: 300 DEG 4MIN/40 DEG/MIN/FINAL 180 DEG

RANGE: G 1, 485 LABEL: N 0, 4.0 QUAN: A 0, 1.0 J 0 BASE: U 20, 3

350512.


100.07

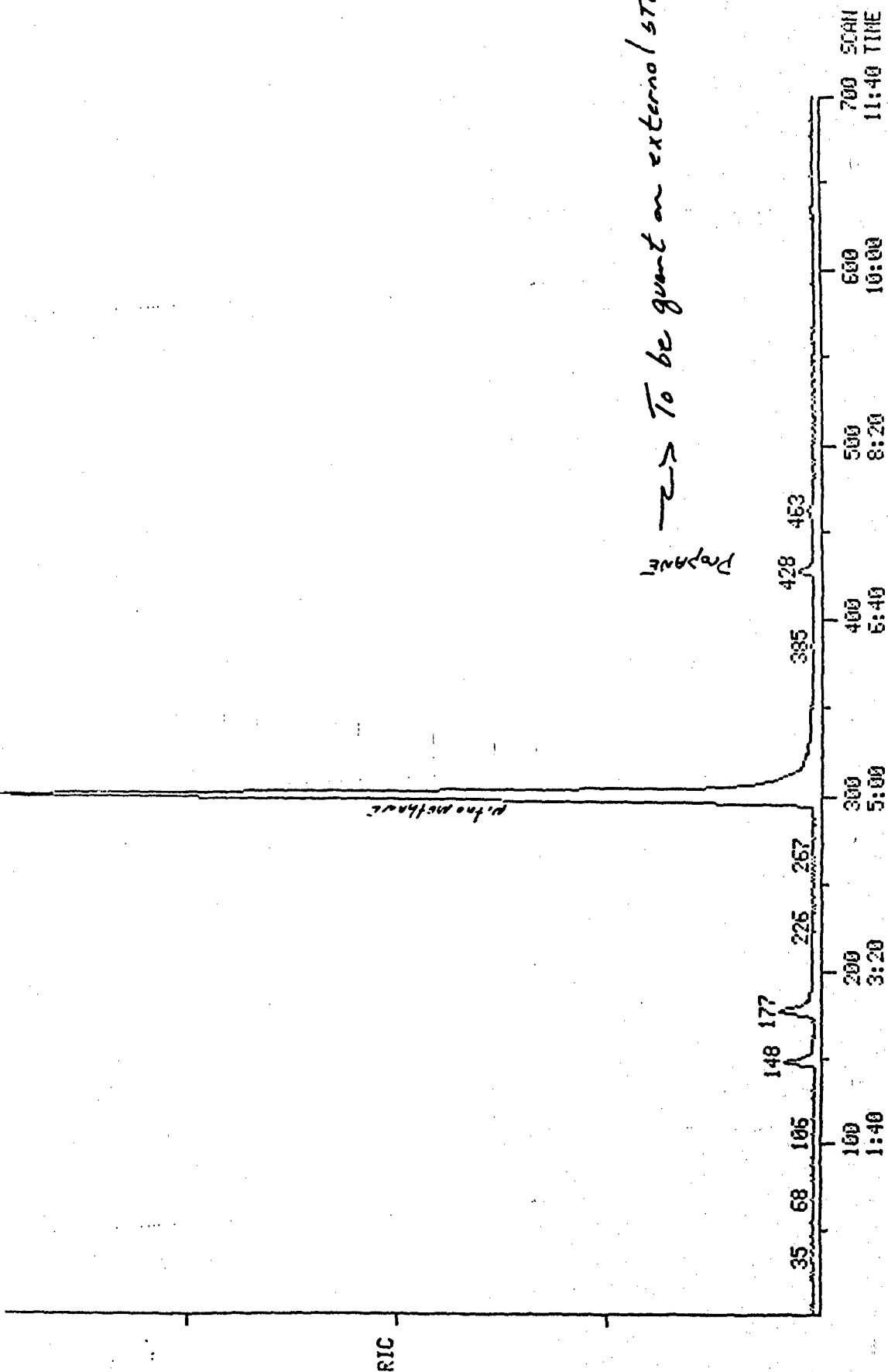

Prepave
→ To be quant on external STD


LIBRARY SEARCH
 06/21/90 22:21:00 + 5:01
 SAMPLE: B.H.:A1-B (SW):1UL(50UL SAMPLE/10ML MECH)
 COND.: 30DEG 4MIN/40DEG/MIN/FINAL 130 DEG
 # 300 TO # 302 SUMMED - # 309 TO # 315 - # 230 TO # 290 X1.00

DATA: 00621001 # 301
 CALI: 006150F1 # 3
 BASE M/Z: 61
 RIC: 2209540.

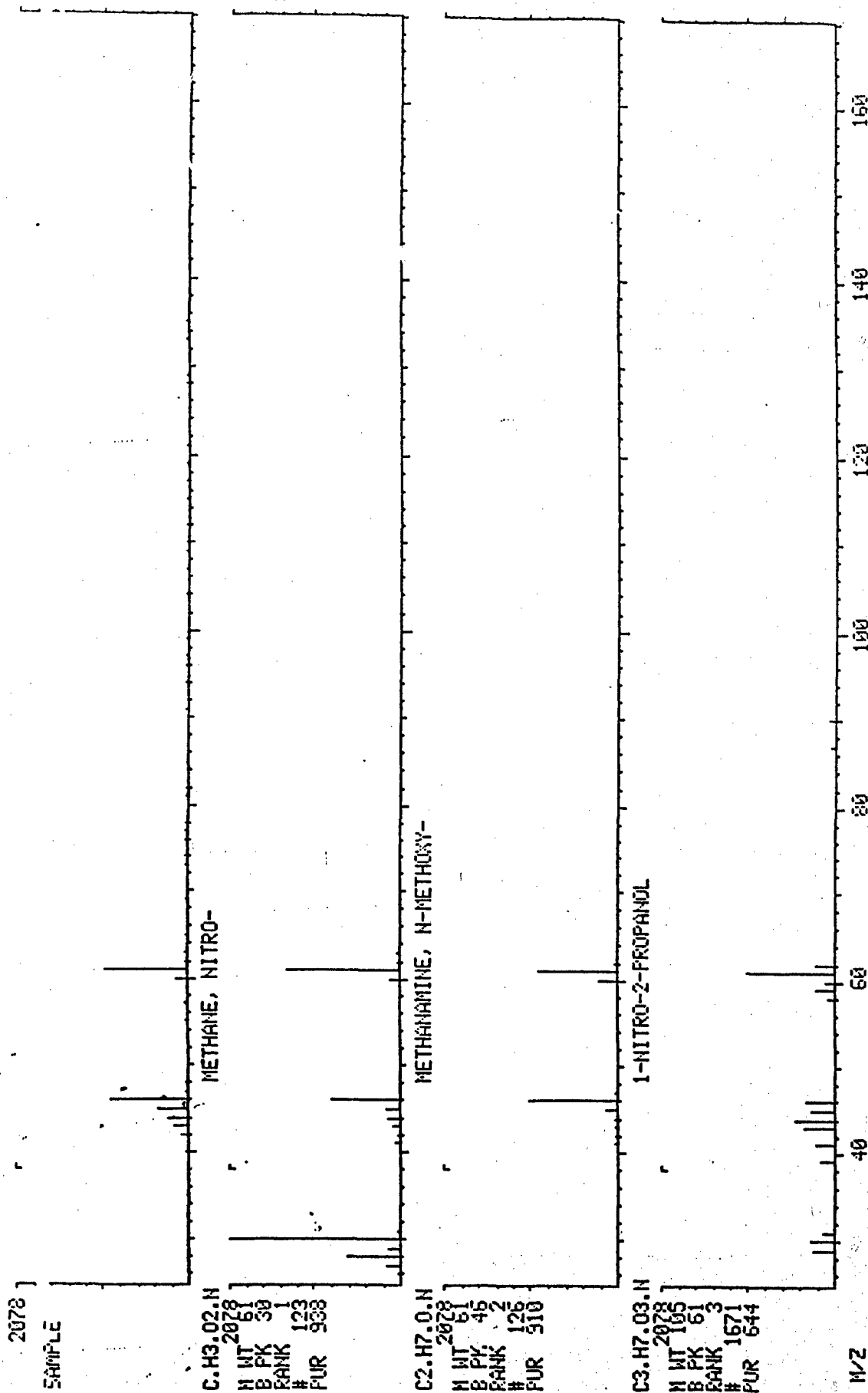

RIC
 05/22/90 0:22:00
 SAMPLE: B.H.: 1UL (100UL/2ML MEDH) 50UL/10ML MEDH
 CONDOS.: 300DEG 4MIN/40DEG/MIN/FINAL 180 DEG
 RANGE: G 1.629 LABEL: N 0.4.0 QUAN: A 0.1.0 J 0 BASE: U 20. 3
 DATA: 00621002 #434
 CALI: 006150F1 #3
 SCHE 1 TO 700

156160.

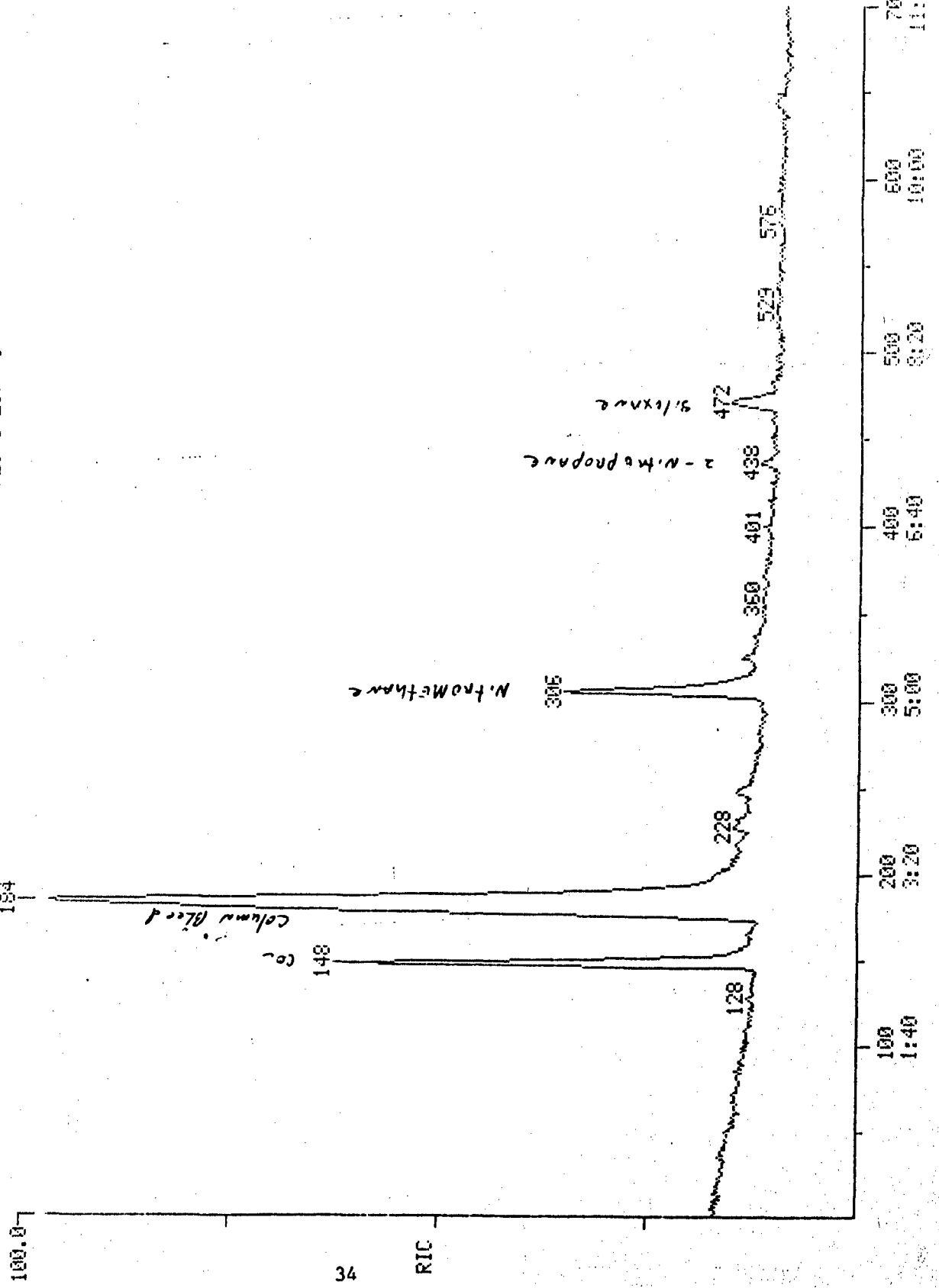

LIBRARY SEARCH

06/22/90 0:22:00 + 5:06

SAMPLE: B.H.: 1UL (100UL/2ML MECH) 50UL/10ML MECH

CONDOS.: 30000 4MIN/40000 MIN/FINAL 100 DEG

304 TO # 303 SUMMED - # 275 TO # 293 - # 333 TO # 353 X1.00

DATA: 00611802 # 305
CALI: 00611801 # 3

DATE: 06/22/90
TIME: 14:25:03

1744

SAMPLE

C.H3.O2.N

1744
N WT 61
B PK 30
RANK 1
123
PUR 908

METRANE, NITRO-

35

C2.H7.O.N

1744
N WT 61
B PK 45
RANK 2
126
PUR 853

METHANAMINE, N-METHOXY-

C3.H7.O3.N

1744
N WT 105
B PK 61
RANK 3
1671
PUR 802

1-NITRO-2-PROPANOL

M/Z

LIBRARY SEARCH


06/22/90 0:22:00 + 7:16

SAMPLE: B.H.: 1UL (100UL/2ML MECH)50UL/10ML MECH

CONDS.: 30DEG 4MIN/40DEG/MIN/FINAL 180 DEG

434 TO # 439 SUMMED - # 448 TO # 456 - # 414 TO # 424 X1.00

DATA: C0621002 # 406
CALI: C0613051 # 3
BASE N 3: 43
R10: 11551.

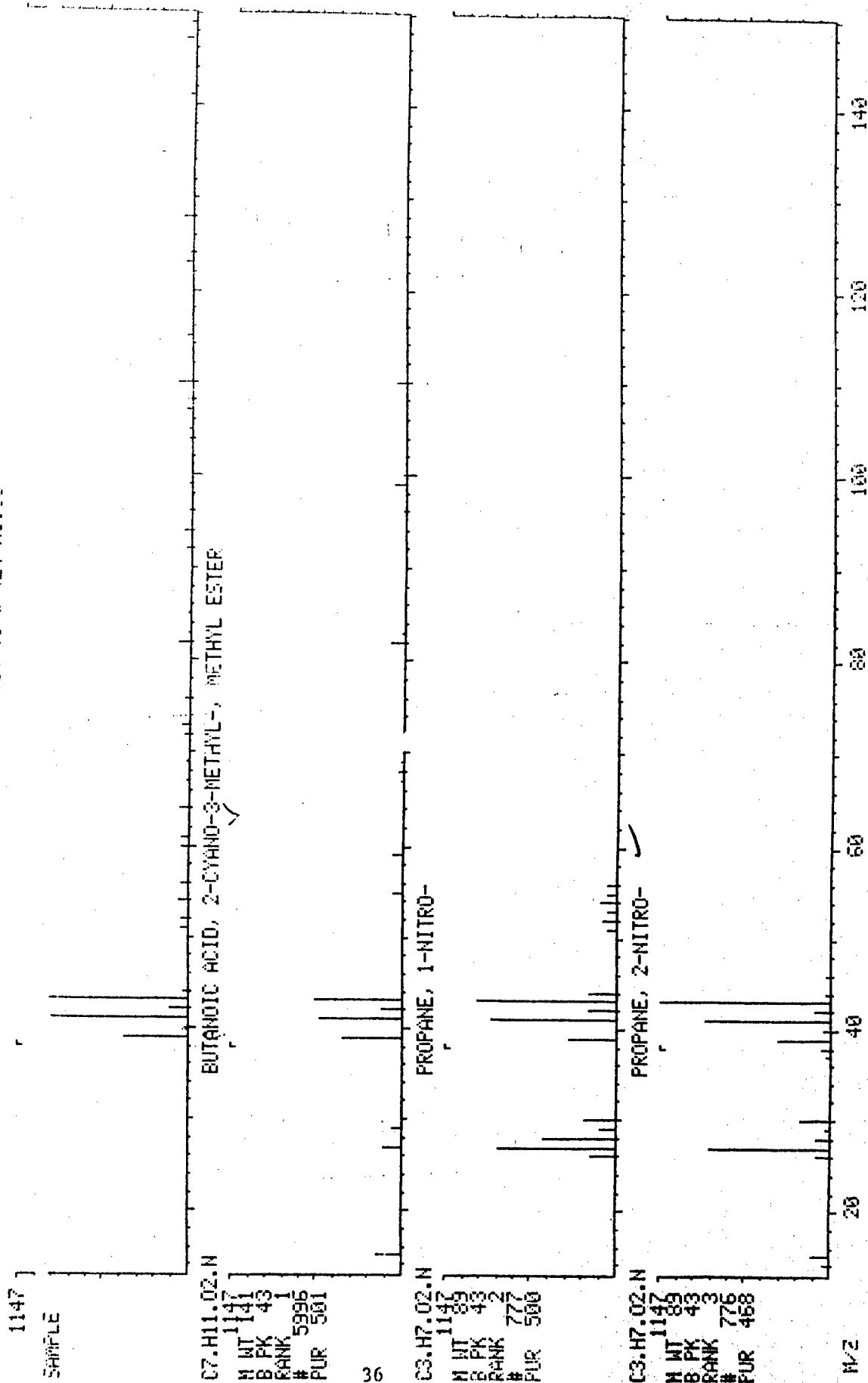

LIBRARY SEARCH

06/22/90 0:22:00 + 7:50

SAMPLE: B.H. :10UL (100UL/2ML MECH)50UL/10ML MECH

COND5.: 30DEG 4MIN/4DEG/MIN/FINAL 100 DEG

468 TO # 473 SUMMED - # 485 TO # 499 - # 434 TO # 448 X1.00

DATA: 00021002 # 473
CALL: 000150F1 # 3

PAGE 171 107
PIC: 41571.

1107

SAMPLE

05.H18.03.S13

CYCLOTRISILOXANE, HEXAMETHYL-

N MT 1107
B PK 207
RANK 1
19998
PUR 893

37

03.H27.03.A5.S13

ARSENIOUS ACID, TRIS(TRIMETHYLSILYL) ESTER

N MT 1107
B PK 207
RANK 2
32875
PUR 763

08.H22.03.S12

DISILOXANE, 1,3-DIETHOXY-1,1,3,3-TETRAMETHYL-

N MT 1107
B PK 207
RANK 3
20022
PUR 726

M/Z

50

100


150

200

250


300

RIC
 05/22/90 0:22:00
 SAMPLE: B.H.: 1UL (100UL/2ML MECH) 50UL/10ML MECH
 COND.: 300 DEG 4MIN/40 DEG/MIN/FINAL 180 DEG
 RANGE: G 1.629 LABEL: N 0, 1.0 J 0 BASE: U 20, 3
 DATA: 00521002 #434
 CALI: 005150F1 #3
 SCANS 1 TO 700
 155150.

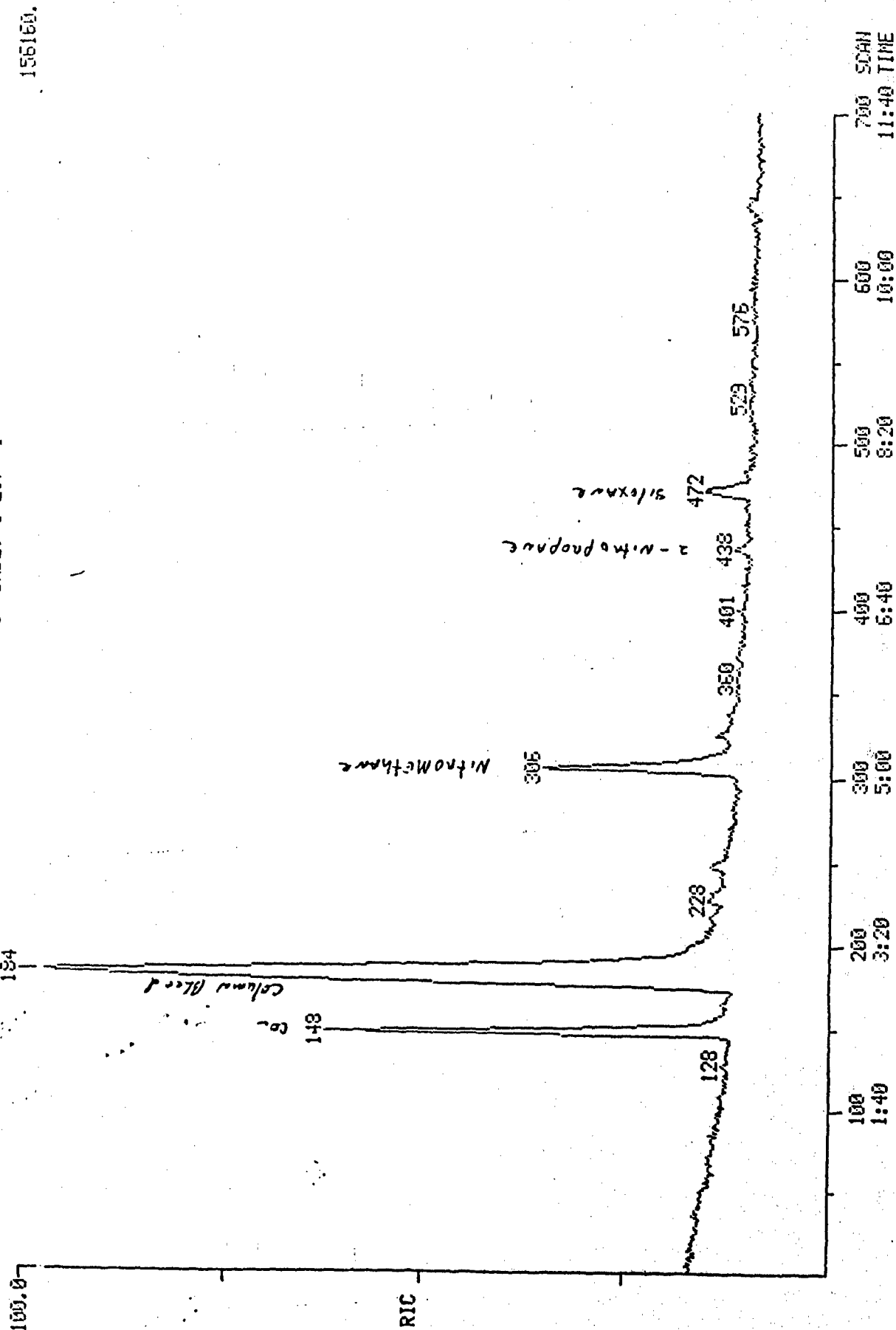

DATA: 00521002 # 305
CALI: 005150F1 # 3
BASE M/Z: 51
RIC: 149503.

LIBRARY SEARCH
05/22/90 0:22:00 + 5:06
SAMPLE: B.H.: 1UL (100UL/2ML MECH) 500UL/10ML MECH
CONDS.: 30DEG 4MIN/4DEG/MIN/FINAL 180 DEG
304 TO # 308 SUMMED - # 275 TO # 293 - # 333 TO # 353 X1.00

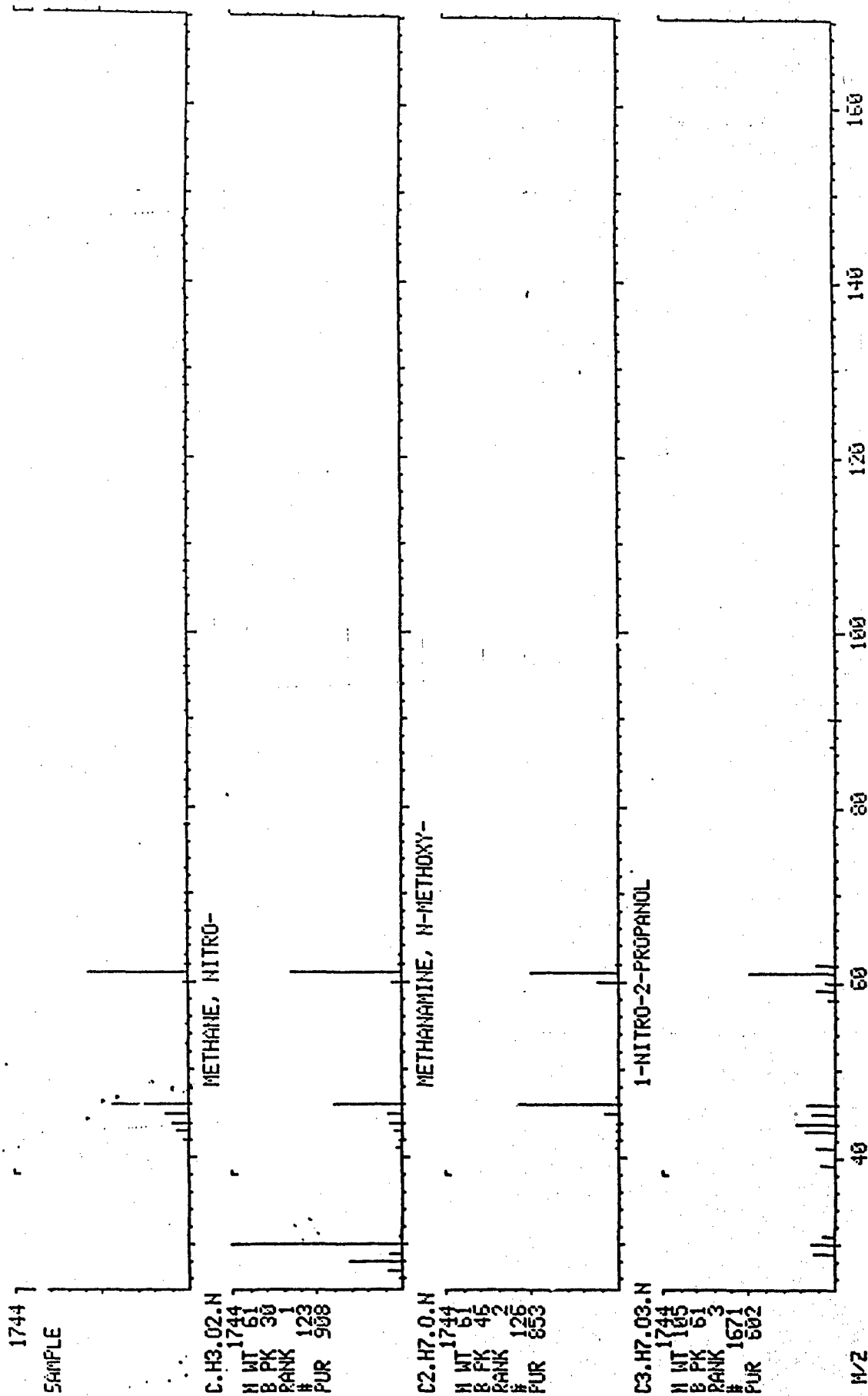

LIBRARY SEARCH

06/22/90 0:22:00 + 7:16


SAMPLE: B.H.: 11UL (100UL/2ML MECH) 500UL/10ML MECH

COND5.: 30DEG 4MIN/40DEG/MIN/FINAL 180 DEG

434 TO # 439 SUMMED - # 448 TO # 458 - # 414 TO # 424 X1.00

DATA: 00521002 # 435
CALI: 005150F1 # 3

EASE M/2: 43
R1C: 11951.

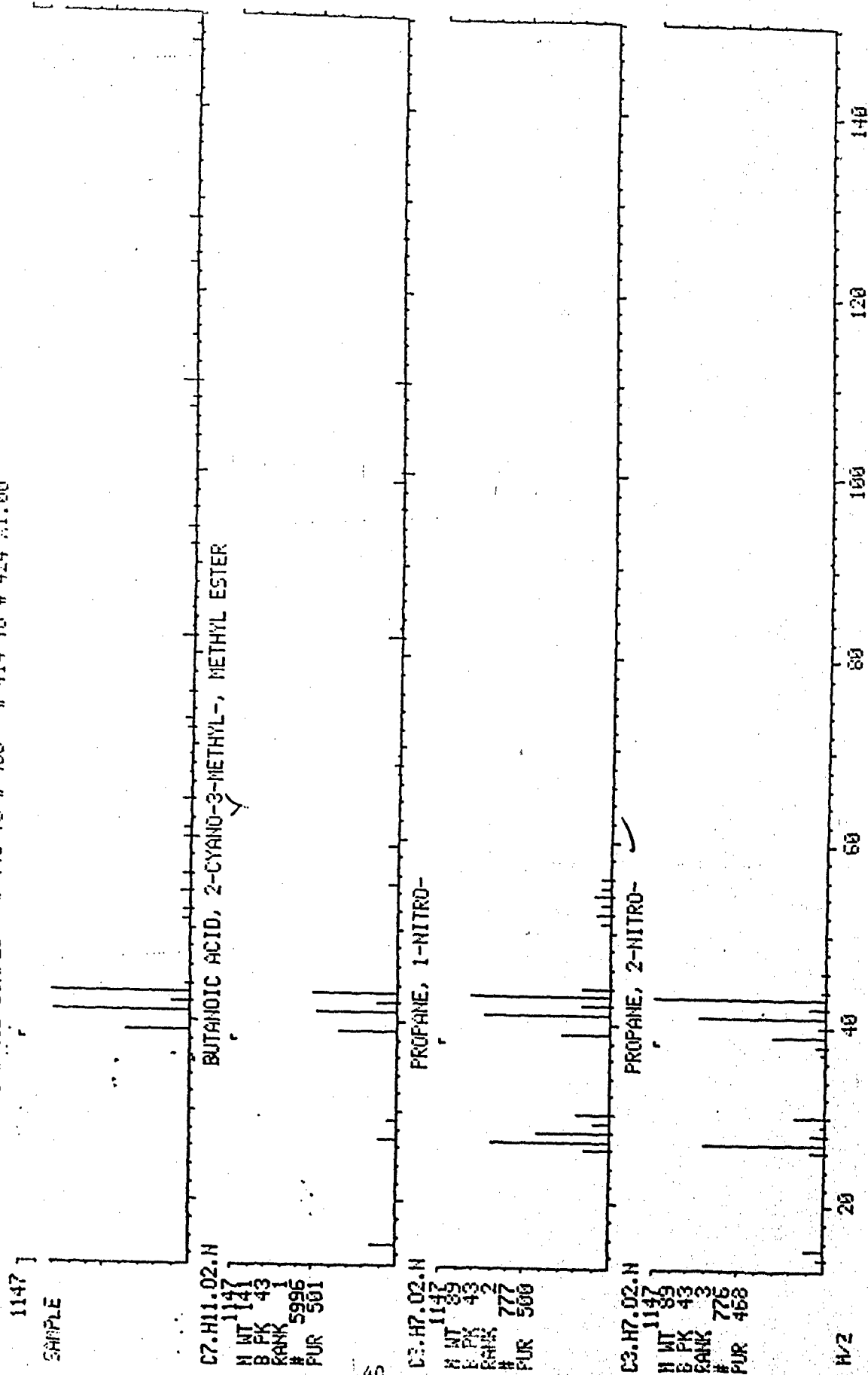

LIBRARY SEARCH

06/22/90 0:22:00 + 7:50

SAMPLE: B.H.:1UL (100UL/2ML MECH)50UL/10ML MECH

COND.: 300DEG 4MIN/4DEG/MIN/FINAL 180 DEG

468 TO # 473 SUMMED - # 485 TO # 499 - # 434 TO # 448 X1.00

DATA: 00521002 # 470
CALL: 005150F1 # 3

BASE M/Z: 287
R1C: 41471.

1107

SAMPLE

C6.H18.03.S13

1107

8 PK 207

RANK 1

19998

PUR 853

CYCLOTRILOXANE, HEXAMETHYL-

C9.H27.03.A5.S13

1107

8 PK 207

RANK 2

32875

PUR 763

ARSENIOUS ACID, TRIS(TRIMETHYLSILYL) ESTER

C8.H22.03.S12

1107

8 PK 207

RANK 3

20022

PUR 726

DISILOXANE, 1,3-DIETHOXY-1,1,3,3-TETRAMETHYL-

M/Z

50

100

150

200

250

300

DATA: 00621052 #1

SCANS 1 TO 700


CALL: 006150F1 #3

SAMPLE: 50NG NITROMETHANE/2-NITROPROPANE STD: 1UL PURGE


COND.: 30DEG 4MIN/40DEG/MIN/FINAL 180 DEG

RANGE: G 1, 492 LABEL: N 0, 4.0 QUAN: A 0, 1.0 J 0 BASE: U 20, 3

105553.

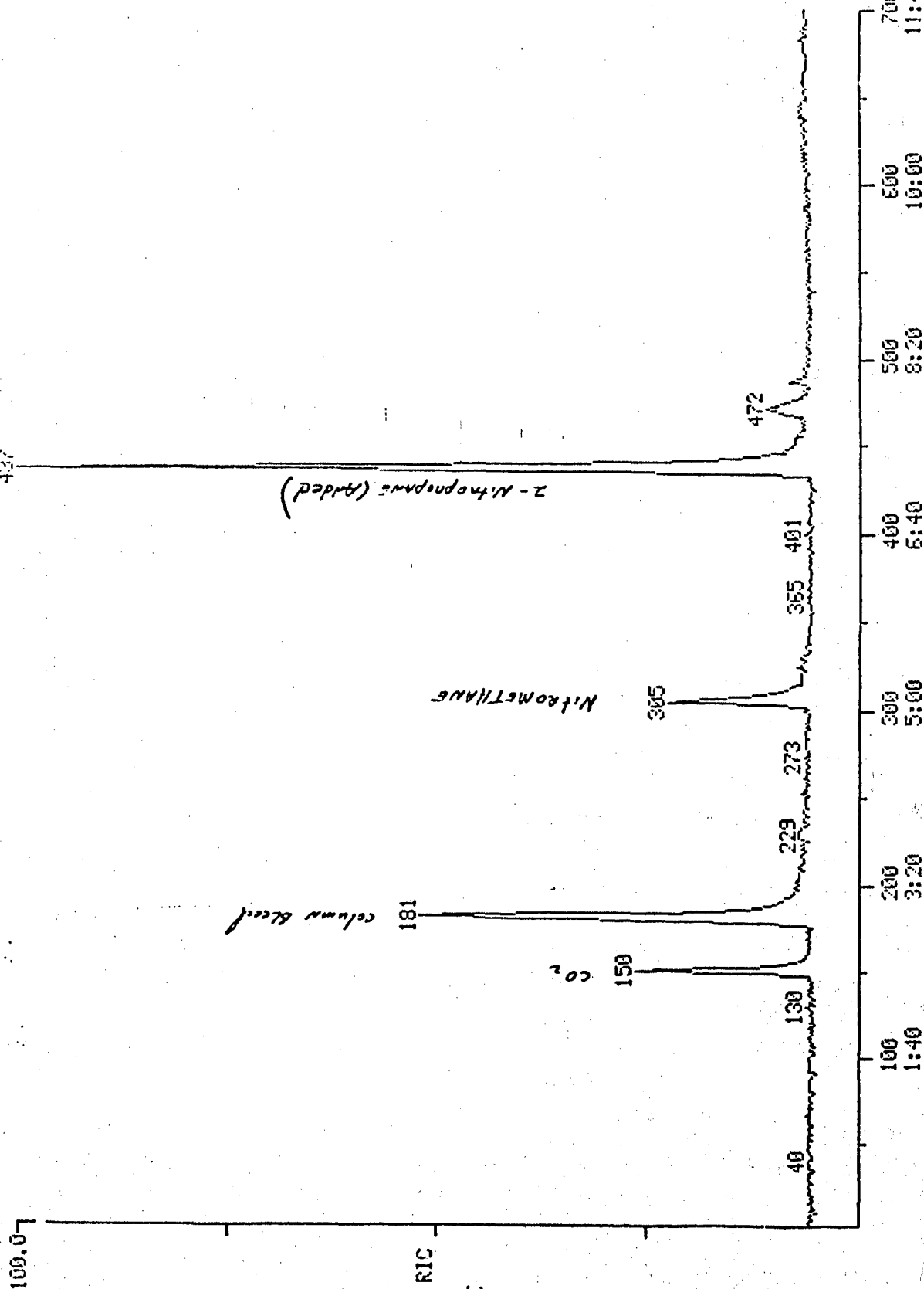

LIBRARY SEARCH
 05/21/90 20:45:00 + 5:06
 DATA: 00521052 # 306
 CALI: 006150F1 # 3
 BASE M/Z: 61
 RID: 53095.
 SAMPLE: 50NG NITROMETHANE/2-NITROPROPANE STD: 1UL PURGE
 COND.: 300DEG 4MIN/40DEG/MIN/FINAL 180 DEG
 # 304 TO # 308 SUMMED - # 315 TO # 321 - # 283 TO # 294 X1.00

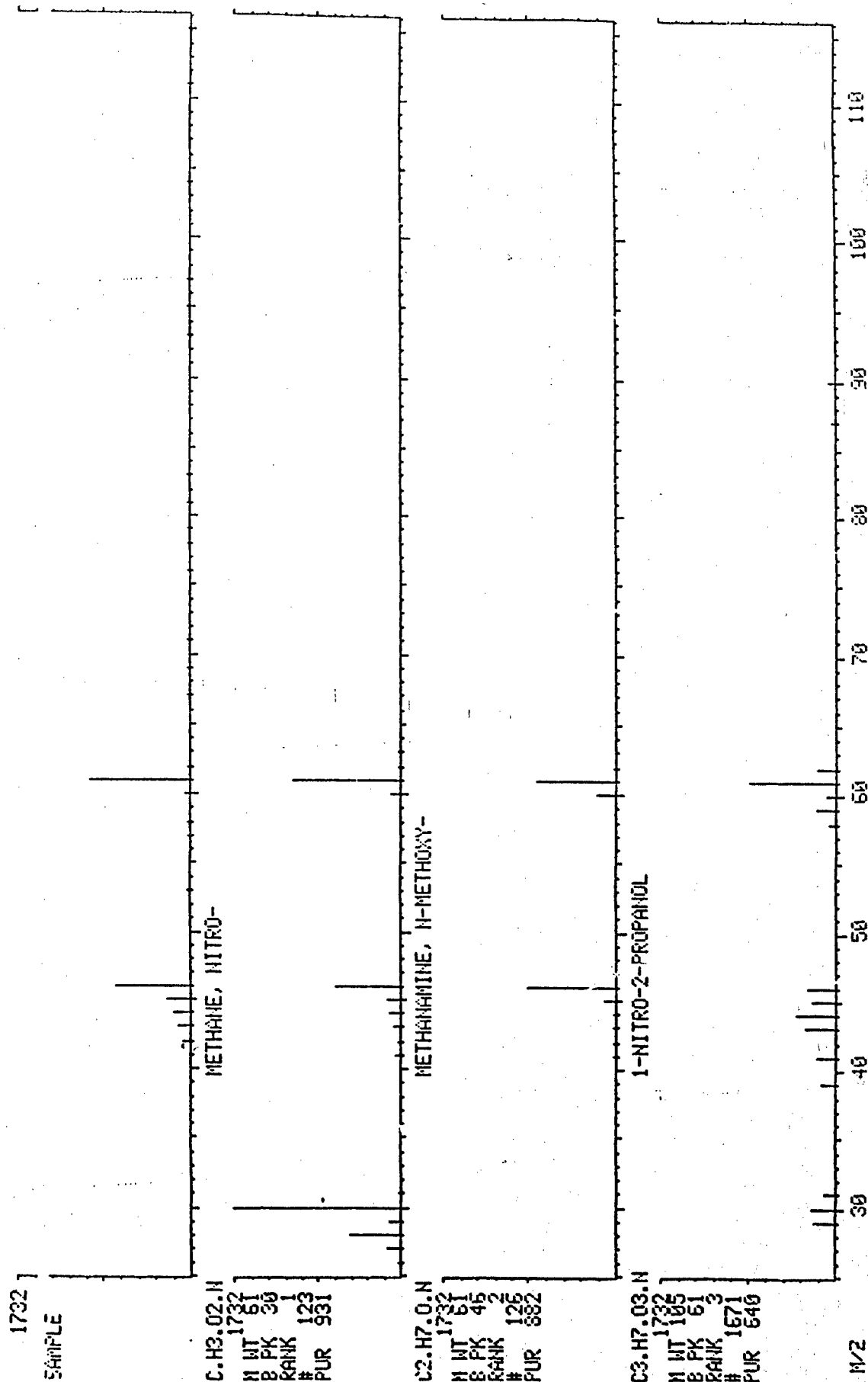

LIBRARY SEARCH

05/21/90 20:45:00 + 7:18


SAMPLE: 50MG NITROMETHANE/2-NITROPROPANE STD: 1UL PURGE

COND.: 300EG 4MIN/40EG/MIN/FINAL 180 DEG

437 TO # 439 SUMMED - # 442 TO # 451 - # 425 TO # 435 X1.00

DATA: 00621052 # 435
CALI: 00615051 # 3

DATE: 43
RIC: 222453.

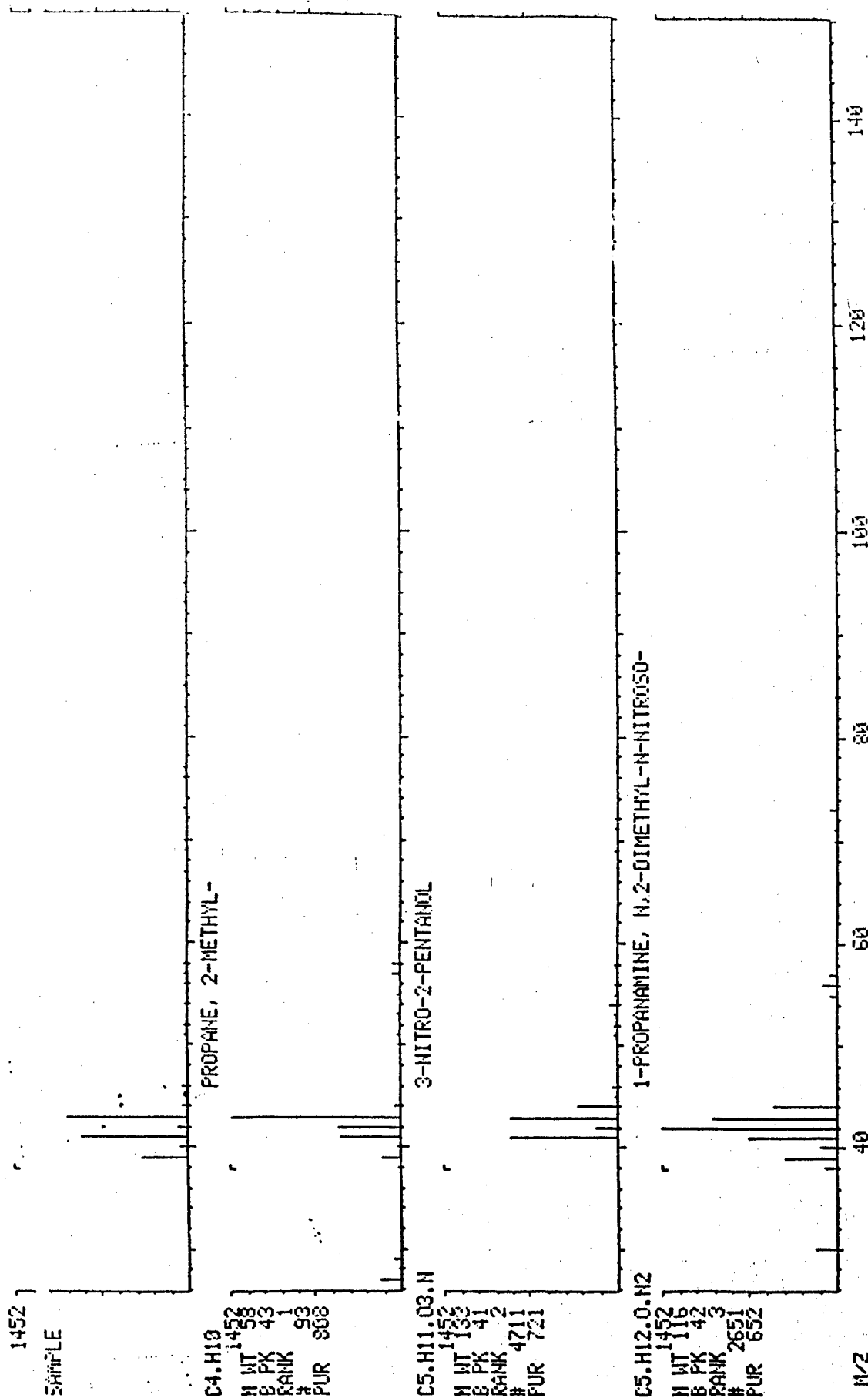

LIBRARY SEARCH
 05/21/90 20:45:00 + 7:17
 SAMPLE: 50NG NITROMETHANE/2-NITROPROPANE STD: 1UL PURGE
 CONDS.: 30DEG 4MIN/4DEG/MIN/FINAL 180 DEG
 # 435 TO # 439 SUMMED - # 443 TO # 449 - # 454 TO # 459 W1.00

DATA: 00621852 # 437
 CALI: 006150F1 # 3

BASE M/2: 43
 RIC: 037407.

